

1 A review of the determination of persistent
2 organic pollutants for environmental
3 forensics investigations

4 **ACA-16-17-Rev**

5 David Megson^{1,2}, Eric J. Reiner^{1,3}, Karl J. Jobst^{3,4}, Frank L Dorman⁵, Mathew Robson ⁶, Jean-François
6 Focant⁷

7 1. University of Toronto, Department of Chemistry, Toronto, ON, Canada,

8 2. School of Science and the Environment, Manchester Metropolitan University, Manchester, UK

9 3. Ontario Ministry of the Environment and Climate Change, Laboratory Services Branch, Toronto,
10 ON Canada M9P3V6

11 4. Department of Chemistry and Chemical Biology, McMaster University, Hamilton, ON, Canada,
12 L8S 4M1

13 5. Penn State University, Biochemistry and Molecular Biology. Forensic Science Program, University
14 Park, PA 16802

15 6. Department of Chemistry, Brock University, St Catharine's, ON, Canada

16 7. University of Liège Organic and Biological Analytical Chemistry , B6c B-4000 Liège Belgium

17

18 **Abstract**

1 The field of environmental forensics emerged in the 1980s as a consequence of legislative frameworks
2 enacted to enable parties, either states or individuals, to seek compensation with regard to contamination or
3 injury due to damage to the environment. This legal environment requires stringent record keeping and
4 defensible data therefore analysis can sometimes be confined to data to be obtained from certified
5 laboratories using a standard accredited analytical method. Many of these methods were developed to target
6 specific compounds for risk assessment purposes and not for environmental forensics applications such as
7 source identification or age dating which often require larger data sets. The determination of persistent
8 organic pollutants (POPs) for environmental forensic applications requires methods that are selective but
9 also cover a wide range of target analytes which can be identified and quantified without bias. POPs are
10 used in a wide variety of applications such as flame retardants, fire suppressants, heat transfer agents,
11 surfactants and pesticides mainly because of their chemical inertness and stability. They also include
12 compounds such as dioxins that can be unintentionally produced from industrial activities. POPs are
13 persistent in the environment, bioaccumulative and/or toxic and therefore require analytical methods that
14 are sensitive enough to meet the low detection limits needed for the protection of the environment and
15 human health. A variety of techniques, procedures and instruments can be used which are well suited for
16 different scenarios. Optimised methods are important to ensure that analytes are quantitatively extracted,
17 matrix coextractables and interferences are removed and instruments are used most effectively and
18 efficiently. This can require deviation from standard methods which can open the data up to further scrutiny
19 in the courtroom. However, when argued effectively and strict QA/QC procedures are followed the
20 development and optimisation of methods based on investigation specific scenarios has the potential to
21 generate better quality and more useful data.

22 **Keywords**

23 Persistent organic pollutants; Environmental forensics; Analysis; Environment

24

1 Introduction

2 The Stockholm Convention on persistent organic pollutants (POPs) targets 26 compounds or compound
 3 groups including the original dirty dozen which comprises of nine organochlorine pesticides,
 4 polychlorinated biphenyls (PCBs) polychlorinated dibenzo-p-dioxins (PCDDs) and polychlorinated
 5 dibenzofurans (PCDFs) [1]. An additional nine pesticides, flame retardants and surfactants were added in
 6 2009, endosulfan was added in 2011, hexabromocyclododecane added in 2013 and hexachlorobutadiene,
 7 pentachlorophenol (and its salts and esters) and polychlorinated naphthalenes in 2015 [1]. There are
 8 currently four compounds/groups: decabromodiphenyl ether (BDE 209), dicofol, short-chain chlorinated
 9 paraffins (SCCPs), and pentadecafluorooctanoic acid (PFOA) its salts and related compounds, that are
 10 currently under review [1] (Table 1).

11 **Table 1 Stockholm Convention POPs List**

Annex	A - Elimination	B- Restriction	C- Unintentional Production
Original List 2001	Chlordane, Dieldrin, Endrin, Heptachlor, Hexachlorobenzene, Mirex, Toxaphene, Polychlorinated biphenyls (PCBs)	DDT	Polychlorinated dibenzo-p-dioxins and polychlorinated dibenzofurans, Polychlorinated biphenyls (PCBs), Hexachlorobenzene
Added May 2009	α -Hexachlorocyclohexane, β -Hexachlorocyclohexane, Chlordecone, Hexabromobiphenyl, Hexabromodiphenyl ether and heptabromodiphenylether, Lindane (gamma-hexachlorocyclohexane) Pentachlorobenzene Tetrabromodiphenyl ether and pentabromodiphenyl ether	Perfluorooctanesulfonic acid (PFOS), and its salts perfluorooctanesulfonyl fluoride (PFOSF)	Pentachlorobenzene
Added May 2011	Endosulfan		

Added May 2013	Hexabromocyclododecane (HBCD)		
Added May 2015	Hexachlorobutadiene Pentachlorophenol and its salts and esters Polychlorinated naphthalenes		Polychlorinated naphthalenes
Under review	Decabromodiphenyl ether (BDE 209), Dicofol, Short-chain chlorinated paraffins (SCCPs), Pentadecafluorooctanoic acid (PFOA), its salts and related compounds		

1

2 This is only a list of priority pollutants and there are over 100,000 chemicals are currently in use or present

3 in consumer products with over 30,000 considered to be in wide commercial use (in consumer products at

4 >907kg/year) [2]. Many of these compounds are persistent, toxic and bioaccumulative and are detected in

5 all types of environmental matrices. In environmental forensics investigations it can be prudent to adopt

6 non targeted methods which can screen for the presence of a large number of contaminants rather than

7 focusing solely on traditional accredited methods that target only a limited number of analytes. Forensic

8 investigations usually involve establishing the source of a contamination event, this can occur from a direct

9 chemical spill, fugitive emissions of chemicals from manufacturing plants, transportation and storage of

10 consumer products and leaking or leaching from final products. A variety of post release processes may

11 alter source signatures of POPs including; volatilization and dispersion, biodegradation, uptake in biota,

12 biotransformation and elimination. This results in environmental samples containing very complex patterns

13 of POPs which can be difficult to interpret. It is therefore important to consider the sample collection,

14 preparation and analysis carefully so that the accuracy and uncertainty of the technique is acceptable to

15 meet the required data quality objectives and the analytical method used is fit for the purpose for which it

16 was intended [3]. This manuscript provides an overview of the available sample preparation methods and

17 analytical techniques. It discusses pertinent issues surrounding method development and identifies points

18 to consider when undertaking an environmental forensics investigation involving POPs.

1 2 Sample preparation

2 2.1 Collection

3 A well designed sampling program is integral to any environmental forensics investigation. There are many
4 factors that can potentially bias data and lead to erroneous results and these need to be understood if they
5 are to be defended in a courtroom. Detailed information on sample collection for environmental forensics
6 investigations has been produced by Mudge [4], Morrison and Murphy [5] and Murphy and Morrison [6].
7 With any investigation it is essential to understand the environmental chemistry of compound in question
8 so that the correct matrix can be targeted. Consideration should also be given to designing an appropriate
9 sample strategy based on the question posed. Commonly asked questions include; what are the potential
10 sources of contamination, what is the extent of the contamination, was the wildlife affected by the
11 contamination or was a regulatory threshold limit exceeded? Each of these questions may involve the
12 collection of different samples for a specific purpose.

13 Numerous samples are gathered during the course of an environmental forensics investigation. These might
14 be collected from a variety of different matrices such as; air particles, gas phases, vegetation, soils, water
15 and animals. In all cases, samples should be clearly labelled and stored in containers appropriate for the
16 analysis being undertaken. All samples should be transported to the laboratory as soon as possible
17 (preferably within 24 hrs) with appropriate chain of custody documentation. Failure to follow appropriate
18 best practice guidance can lead to errors in the data and may ultimately result in the findings being dismissed
19 from the courtroom. It is clearly preferable to follow the best practice guidance, however failure to follow
20 best practice guidance does not necessarily mean that the data should be completely disregarded [7]. There
21 have been cases in the U.S. such as *People v. Hale* (56Cal.App3d [1994]) where deviations from the best
22 practice guidance did not preclude the introduction of the information as evidence [8].

1 2.2 Extraction

2 Extraction of the sample is a critical step that is required to separate the compounds of interest from the
3 bulk matrix. Prior to extraction, internal standards are typically added to the sample to allow for
4 quantification, quality control (QC) and to measure the extraction efficiency and recovery rate. Isotopically
5 labelled analogues of the target compounds are preferred as they behave the same way as the native
6 compounds in the sample which results in the most accurate and precise data for calibration and
7 quantification [9]. However, isotope labelled solutions can be very expensive and in some cases there may
8 not be labelled standards available. In these cases a representative compound with similar properties to the
9 target compound can be used to calculate recovery rates. The accuracy and precision of the method can be
10 strongly dependent on the number of internal standards and type of calibration used. In environmental
11 forensics investigations samples can be obtained from a variety of different matrices, each of these may
12 require a different extraction procedure to efficiently separate the compounds of interest.

13 2.2.1 Air, aqueous and liquid matrices

14 Aqueous and other liquid samples like biological fluids have classically been extracted using liquid/liquid
15 extraction. If samples contain particles, they can be filtered to isolate the precipitates/particles and the two
16 fractions either determined separately or recombined at a later stage depending on the aims of the project.
17 This is particularly important for hydrophobic compounds such as BDE 209 that may be present almost
18 exclusively in the particulate phase. An alternative to liquid/liquid extraction that has become more
19 common is solid phase extraction (SPE). This technique uses a stationary phase / resin in an extraction
20 cartridge or disc to extract non-polar analytes like dioxins from polar liquids. Extraction of samples by SPE
21 [10-12] also allows for the particles to be trapped on top of the extraction disc or bed. Quantitative elution
22 of the particles and disk can be done in several steps to separate several compound classes, or in one step
23 for many analytes which can significantly reduce solvent usage and analytical time [9]. There are several
24 other methods that have been successfully used for the extraction of POPs such as; stir-bar sorptive
25 extraction (SBSE) with thermal desorption directly into the gas chromatograph [13], hollow fiber liquid-

1 phase microextraction (HF-LPME) [14], Solid phase micro extraction (SPME) [15], and passive samplers,
2 including semipermeable membrane devices (SPMD) or polyethylene strips [16]. Air sampling for POPs is
3 traditionally conducted using polyurethane foam (PUF) or XAD resin [17], or a mixture of both to capture
4 the more volatile gas phase analytes, combined with a Teflon or glass fibre filter (GFF) to capture those
5 analytes found predominantly in the particulate phase. More recently researchers have also used passive
6 samplers such as low-density polyethylene passive samplers (LDPE-PAS) to sample both the aquatic and
7 atmospheric environment [16]. PAS have the benefit of little to no infrastructure requirements, easy
8 extraction, potentially long deployment times and high capacities for POPs. PAS have been used in a
9 number of studies for the source tracking of PCBs [18-20] and OCPs [21, 22].

10 2.2.2 Solid matrices

11 Solid matrices (e.g., soil/sediment, biota or vegetation) have classically been extracted using Soxhlet
12 extraction or an automated version of Soxhlet called Soxtec. Pressurized liquid extraction (PLE – also
13 known as accelerated solvent extraction (ASE)) [23-25], microwave assisted extraction (MAE) [26, 27] and
14 sonication have also be used. PLE and MAE are automated techniques that subject the samples to elevated
15 temperatures and pressures that can result in faster more efficient extractions than Soxhlet and sonication.
16 The use of PLE for the extraction of POPs in environmental samples has recently been reviewed by
17 Vazquez-Roig and Picó [28]. Supercritical fluid extraction (SFE) [29, 30] is a procedure that uses
18 supercritical carbon dioxide and requires no solvent for extraction and only a very small amount of solvent
19 to trap the analytes or elute them from the carbon or C18 trapping material. Solid matrices can also be
20 digested by using an acid or base and then extracted using SPE or liquid/liquid extraction [31, 32], however
21 care must be taken when using an acid or base as they may breakdown potential compounds of interest such
22 as BEHTBP and other susceptible HFRs. Many of the techniques discussed in section 2.2 are relatively
23 expensive as they require a lot of time and need to be conducted by an experienced scientist to produce
24 acceptable results. QuEChERS (quick, easy, cheap, effective, rugged, safe) provides a cheaper alternative
25 where the cost per sample can be less than \$10. Many specific QuEChERS methods have been developed

1 for different applications based on the original method of Anastassiades et al. [33]. The sample is extracted
 2 in the presence of a dispersant in a disposable tube providing a relatively clean injection ready extract,
 3 although additional cleanup steps are often used. Although initially developed for the determination of
 4 veterinary drugs the method has modified to determine non-polar analytes like PAH and PCBs and is an
 5 excellent choice for initial sample screening [34-36]. The use of QuEChERS for the determination of
 6 persistent organic pollutants in environmental matrices has been reviewed [37, 38].

7 A summary of available extraction techniques for the different matrices is presented in Table 2. It provides
 8 a brief overview of each technique and lists a couple of pertinent references which can be referred to for
 9 further information.

10 **Table 2. Summary of extraction techniques**

Method	Main Matrix(es)	Brief overview	Example references
Solid Phase Extraction (SPE)	Aqueous / Liquid	A sample is passed through a stationary phase / resin in an extraction cartridge or disc which extracts analytes from liquids. These can then be collected with a solvent.	[9,10,11,12]
Stir-bar sorptive extraction (SBSE)	Aqueous / Liquid	A small bar coated with sorbent is added to the liquid mixture and stirred to extract analytes. These can then be collected with a solvent.	[13]
Hollow fibre liquid phase microextraction (HF-LPME)	Aqueous / Liquid	Sample is passed into a small hollow fibre coated with adsorbent where it is retained. These can then be collected with a solvent.	[14]
Solid phase micro extraction (SPME)	Aqueous / Liquid	An adsorbent coated fibre is inserted into the sample to extract analytes. These can then be collected with a solvent or the fiber inserted directly to the injector.	[15]
Semipermeable membrane devices (SPMD)	Aqueous / Liquid	The sampler is normally placed in a watercourse and left for a period to adsorb analytes. These can then be extracted from the membrane with a solvent.	[16]
Low density polyethylene strips	Aqueous / Liquid & Air	A strip of polyethylene is placed in a watercourse or outdoor sample case and left for a period to adsorb analytes. These	[16, 19, 21, 22]

(LDPE-PAS or LOPE)		can then be extracted from the membrane with a solvent.	
Polyurethane foam (PUF) and XAD resin	Air	Sample is pumped through a foam / resin disk/filter which absorbs gaseous phase analytes. It can be used in combination with a Teflon or glass fibre filter to trap particulates. These can then be collected with a solvent.	[17,18, 20]
Pressurised liquid extraction (PLE) & Accelerated solvent extraction (ASE)	Solid	Sample is placed in a cartridge and solvent passed through under a designated temperature and pressure to extract analytes.	[23-25, 28]
Microwave assisted extraction (MAE)	Solid	Samples are placed in a temperature controlled vessel and microwaves used to digest sample and extract analytes into a solvent.	[26, 27]
Supercritical fluid extraction (SFE)	Solids	Samples are placed in a temperature controlled vessel and pressurised with supercritical CO ₂ to dissolve the sample, the gas is then depressurised causing the analytes to precipitate in a collector with a small amount of solvent.	[29, 30]
Saponification / acid digestion and extraction	Solid	Sample is digested with an acid or base, the resultant solution can be separated by liquid/liquid extraction or SPE. These can then be collected with a solvent.	[31,32]
QuEChERS (quick, easy, cheap, effective, rugged, safe)	Solids	A solvent (typically acetonitrile) is added to the sample along with water and salts (typically magnesium sulphate and sodium acetate). This is shaken and the salts precipitate out leaving the analytes in the solvent.	[33, 34-38]

1

2 2.3 Sample extract cleanup

3 Most extraction procedures for POPs are quantitative and able to extract a variety of organic compounds as
4 well as other organic matrix co-extractables. Therefore, extract cleanup procedures are required to remove
5 the potentially large quantities (up to a gram or more) of organic material and interfering compounds, but
6 still retain as much of the desired analytes as possible. Most clean up procedures involve combinations of
7 silica, alumina, Florisil® and carbon adsorbents or size exclusion materials (e.g. gel permeation

1 chromatography) to remove matrix and interfering compounds [39-44]. If required more aggressive clean
 2 up techniques can be used such as the use of sulphuric acid to clean up extracts for PCB determination
 3 (EPA Method 3665A). However, care must be taken as the method may breakdown many pH sensitive
 4 analytes including the pesticides Aldrin, Dieldrin, Endrin, Endosulfan (I and II), and Endosulfan sulfate as
 5 well as halogenated flame retardants such as BEHTBP and EHTBB. There are several on-line cleanup
 6 systems that have also been used to provide an automated high throughput alternative [45, 46].

7 In forensics cases it is very important to ensure that the absence of a key compound from a sample is not
 8 due to the sample preparation step. Selecting the appropriate concentration technique is critical as using an
 9 inappropriate procedure could result in significant losses and bias. For further reading Reiner et al. [9],
 10 Lambropoulou [47], Fidalgo-Used et al. [48], Rawa-Adkonis et al. [49], Tang, [50] and Xu et al. [51] have
 11 all reviewed a variety of extraction and clean up techniques for environmental samples.

12 3 Sample analysis

13 3.1 Chromatographic separation

14 In order to achieve accurate quantification, compounds must be fully resolved from each other. In a number
 15 of cases this may be possible by using selective detection. While several chromatographic techniques are
 16 potentially useful, gas chromatography (GC) and high performance liquid chromatography (HPLC) are
 17 clearly the most widely used techniques for environmentally-relevant separations. Since the early 2000's,
 18 ultra-high performance liquid chromatography (UHPLC) has also been successfully used for many
 19 environmental separations, mostly in an effort to increase laboratory throughput [9]. An overview of the
 20 different separation techniques is presented in Table 3.

21 **Table 3. Summary of separation techniques**

Technique	Advantages	Disadvantages	Analytes investigated using the method and selected references
-----------	------------	---------------	--

Liquid chromatography	<ul style="list-style-type: none"> •wide range of potential analytes including non-volatile, thermally labile, polar and ionic analytes •short run times •variety and complementary nature of potential ionisation techniques 	<ul style="list-style-type: none"> •Lower number of theoretical plates gives poorer separation potential compared to GC 	OC Pesticides [57,58, 113, 117, 124] PBDEs [135, 136, 137] HBCD [55, 56, 60-62, 135, 136, 137] Perfluorinated chemicals [64-68, 141, 142]
Gas chromatography	<ul style="list-style-type: none"> •Good separation potential 	<ul style="list-style-type: none"> •Restricted to use on more volatile compounds 	OC Pesticides [45, 57, 112-116, 118, 121-123] PCBs &/or PCNs [45, 121, 123-128] PCDD/F [99, 123, 129-131] PBDEs [120, 127, 128, 132, 134, 135, 136, 138, 139] HBCD [138, 139] Perfluorinated chemicals [51]
Multidimensional chromatography	gas <ul style="list-style-type: none"> •Excellent separation potential 	<ul style="list-style-type: none"> •Restricted to use on more volatile compounds •Large file size and more data processing often required 	OC Pesticides [86] PCBs &/or PCNs [32, 82] PCDD/F [89] PBDEs [133]

1

2 3.1.1 Liquid chromatography

3 In HPLC separations the difference in retention is a complex process involving equilibria of the analytes
 4 between the stationary phase ligand, the mobile phase composition, and possible adsorption forces to the
 5 particle. Since this competitive equilibria is complex, it is not easy to determine how to adjust the chemical
 6 variables in order to improve separation. Modelling programs are available that aid in the optimization of
 7 the chemical variables, however these can be expensive. A cheaper option can be to use predictive models
 8 based on tabulated data. Tabulations based on the hydrophobic subtraction model
 9 (<http://www.usp.org/app/USPNF/columnsDB.html>) are particularly useful in selection and also in cross-
 10 over of column chemistry for method development [9]. Reverse phase liquid chromatography is the most
 11 widely used separation technique for environmental samples. This is usually performed using octadecyl
 12 (C18)-bonded silica and octyl (C8)-bonded silica. Other bonded phase materials on polymeric, monolithic,

1 non-porous and superficially porous supports and special aqueous phase materials have also become
2 available [52]. The use of LC-MS in the determination of contaminants in the environment and foodstuffs
3 has recently been reviewed by Aceña et al. [53] and Hird et al. [54].

4 The key benefits of LC separation in environmental forensic investigations are its wide range of potential
5 analytes (e.g. non-volatile, thermally labile, polar or ionic analytes), short run times (often less than 15
6 mins), and variety and complimentary nature of potential ionisation techniques (e.g. ESI, APCI, APPI). A
7 good example of the benefits of LC analysis in environmental forensic investigations is shown by Zhang et
8 al. [55] who determined the diastereoisomer profiles of HBCD in soils from coastal regions in Northern
9 China and identified different industrial sources of HBCD to the environment. Similarly, Managaki et al.
10 [56] used a HBCD diastereoisomer specific LC-MS/MS method to determine the likely sources of HBCD
11 to the Tsurumi river catchment flowing through the Tokyo metropolitan area.

12 The combination of fast speed of analysis achieved by liquid chromatography in combination with high
13 sensitivity/selectivity with tandem MS detection often makes LC the preferred choice for the determination
14 of pesticides [57, 58] although similarly short (<10 min) run times have recently been achieved by gas
15 chromatography [59]. Most POPs have physiochemical properties that make them well suited to analysis
16 by gas chromatography. The improved degree of separation that is achievable by gas chromatography
17 allows for a more detailed congener specific detection which often makes analysis by gas chromatography
18 preferable in environmental forensics investigations. However, LC is the preferred choice for polar non-
19 volatile and/or thermally labile POPs, such as hexabromocyclododecane (HBCD) [60, 61]. Sales et al. [62]
20 recently proposed the use of gas chromatography with atmospheric pressure chemical ionisation for initial
21 HBCD screening and quantification and compared the results against previously used GC and LC methods.
22 However, there are three main HBCD isomers (α , β , and γ) which interconvert at elevated temperatures
23 during analysis by GC-MS and therefore only total HBCD should be determined this way [63].

24 The perfluorinated chemicals identified by the Stockholm Convention are predominantly determined by
25 liquid chromatography [64]. PFOSF (a major precursor of PFOS) determination by GC has been attempted,

1 although it was not well retained by a variety of GC columns (DB-5, DB-1 and MEGA WAX) [51]. The
2 majority of analytical methods are targeted and involve MS/MS, however more recently methods have been
3 developed to screen for the presence of other per and polyfluorinated alkylated substances (PFASs) [65-
4 68].


5 3.1.2 Gas chromatography

6 When using gas chromatography the separation of compounds in a mixture occurs chiefly due to the
7 different boiling points of each compound and interactions with the stationary phase of the column.
8 Understanding how the stationary phase chemistry effects GC separation is fairly easy to predict and
9 calculations can be performed using commercially available software packages such as “Pro ez-GC”
10 (Restek Corporation, Bellefonte, PA). These programs can be exceptionally useful for both column
11 selection and the optimization of a separation, even for compounds that are not found in the included
12 libraries [9]. By optimising the resolution of compounds as well as possible within the chromatographic
13 space it places less demand on detection systems, and thus improves the overall quality of the total analysis.
14 This can be particularly useful in environmental forensics investigations as the ability to resolve more
15 compounds will produce a more detailed chemical fingerprint which increases the power of statistical tests
16 that are often used for source identification. There are many different chromatographic phases available
17 from a variety of manufacturers. A list of selected available phases and related applications for the
18 determination of POPs are summarised in the supplementary information (SI 1). Selection of a general
19 purpose column (such as a 5% diphenyl dimethylpolysiloxane) can be a good initial step for screening
20 samples to identify potential contaminants of concern or compounds of diagnostic value. This can then
21 inform further phases of analysis which may involve the use of a more species specific column to focus on
22 compound specific analysis (such as a propriety PCB column) or chiral analysis (such as cyclodextrin)
23 which can be used to answer specific environmental forensics questions [32].

1 3.1.3 Multidimensional chromatography

2 Comprehensive two dimensional gas chromatography (GCxGC) was developed at the end of the 20th
3 century by the late John Philips [69]. The technique has several advantages over single dimensional
4 chromatography as it can significantly increase peak capacity (selectivity), increase sensitivity and reduce
5 analysis times by eliminating ejections on multiple chromatographic phases. Detailed reviews on the
6 development and applications of multidimensional chromatography have been produced by Adahchour et
7 al. [70], Bordanjandi et al. [71], Mondello et al. [72], Murray [73], Meinert and Meierhenrich [74], Mostafa
8 et al. [75], Marriott et al. [76], Ramos [77] and Seeley and Seeley [78]. In GCxGC, two GC columns of
9 different phases are connected by the modulator, a thermal or valve controlled device, which traps
10 compounds from the first column and then re-injects them in a tight band onto the second column usually
11 of a complimentary phase where further separation occurs. The best chromatography is achieved when
12 compounds are passed through two columns containing different stationary phases, which results in two
13 degrees of separation based on different physiochemical properties, e.g. boiling point and polarity. Ideally
14 this results in an orthogonal separation where isomers elute in bands at about a 45 degree angle (Figure 1).
15 This also allows several groups of POPs to be identified in the same analytical run and in most cases can
16 provide an extra degree of separation as the target analytes are separated from any remaining sample matrix.
17 One of the greatest challenges of this technique is to be able record a minimal number of measurements
18 across the very narrow peaks produced in the modulation process which are approximately 400 ms wide.
19 In order to accurately define a second dimension GC peak, about 20 spectra per second are required to
20 obtain seven to ten data points. Less than 7 measurements will not be representative of a Gaussian-shaped
21 peak and can therefore result in significant quantitative bias [9].

22

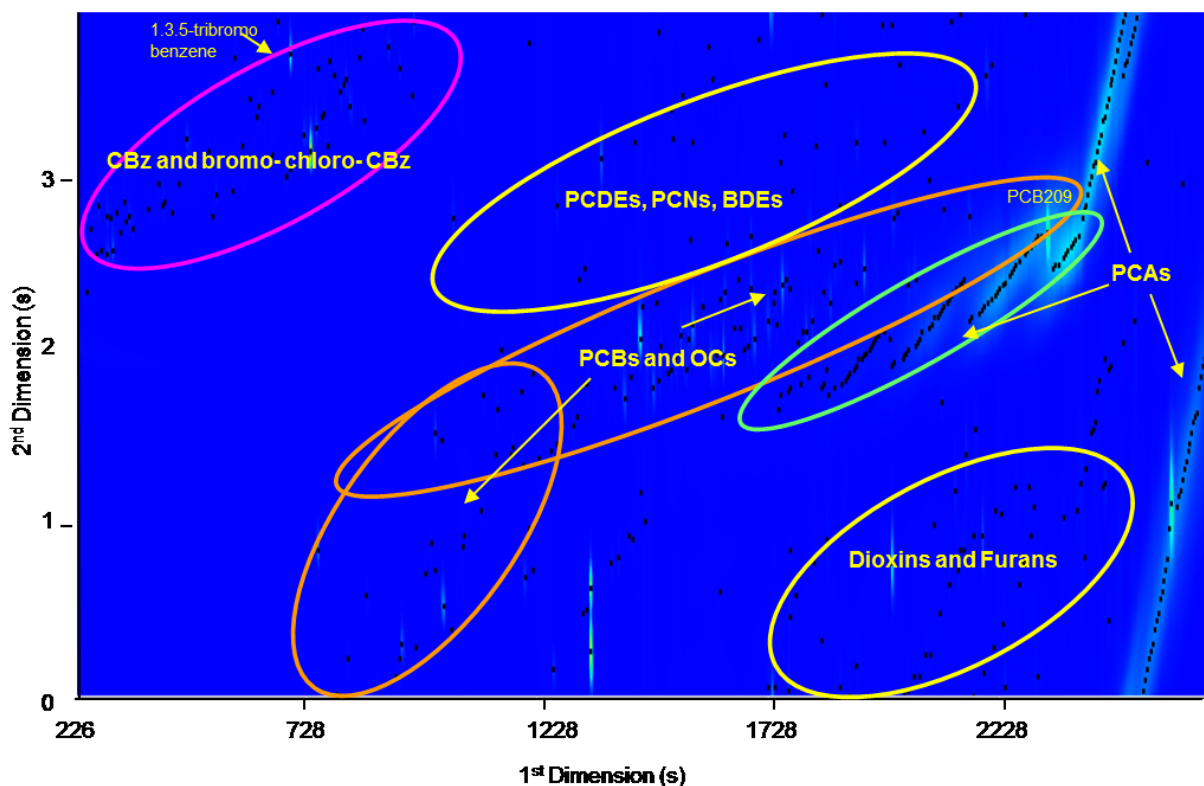

1
 2 **Figure 1.** GCxGC-ECD chromatogram for POPs determination of a sediment sample (reproduced with
 3 permission from Reiner et al. [9]).

4 There are several major advantages to analysis by GCxGC over one dimensional techniques. It can
 5 eliminate the need for fractionation or multicolumn analysis and the modulation process produces much
 6 narrower and taller chromatographic peaks which can increase signal-to-noise ratios (sensitivity) by up to
 7 an order of magnitude. As well as being utilised for targeted analysis GCxGC has been used for the
 8 screening of non-target compounds [79-86] The complex patterns detected in the GCxGC chromatograms
 9 can also be used as fingerprints to characterize sources of contamination. The detailed congener specific
 10 data generated by GCxGC makes it very well suited for environmental forensics investigations involving
 11 complex biological matrices [32, 87]. GCxGC techniques are now at the stage where they can be used for
 12 routine analysis, an example of this is the determination of PCBs, organochlorines and chlorobenzenes by
 13 GCxGC- μ ECD using the Ontario Ministry of the Environment and Climate Change method MOE E3487
 14 [86]. Despite the unparalleled peak capacity; the successful determination of complex environmental,

1 petrochemical, and biological samples; and substantial advancements in instrumentation, the use of GCxGC
2 in the forensic domain has remained relatively scarce since its introduction in 1991 [88]. A review of
3 GCxGC applications in forensics investigations was undertaken by Sampat et al. [88], who recorded only
4 36 peer reviewed articles by 12 research groups/collaborations between 2001 and 2015.

5 3.2 Analyte detection

6 Once the compounds of interest have been separated there are a variety of options available for detection
7 and quantification. Sensitivity, selectivity and speed of analysis are the main considerations when selecting
8 the detector. Nowadays, simple detectors such as flame ionization and electron capture detectors (FID and
9 ECD) are mainly used for quick screening, although these techniques can still prove powerful tools,
10 especially when coupled to GCxGC [86, 89]. ECD detectors are especially sensitive to halogenated
11 compounds and so can provide comparable detection limits to many mass spectrometers for halogenated
12 compounds such as PCBs (Table 4). Recent advances in mass spectrometry have resulted in commercially
13 available instruments that are highly sensitive, selective, automated, and simple to operate with reduced
14 cost and maintenance. As a consequence, bench-top quadrupole and time-of-flight mass spectrometers have
15 in some case supplanted FID and ECD for both screening and quantitative applications. More sophisticated
16 tandem mass spectrometers (MS) and high resolution MS instruments have become more accessible and
17 these can be found in some form in most analytical laboratories. Mass Spectrometers are a versatile detector
18 which offers excellent sensitivity, dynamic range and also offers the key advantage of being able to use
19 mass labelled internal standards, which can significantly increase accuracy and precision [9]. For qualitative
20 work, no other technique can provide the same wealth of structure information of little more than a few
21 picograms of a chemical compound. Detailed reviews on the instrumental determination of POPs have been
22 produced by Guo and Kannan [90], Van Leeuwen and de Boer [40], Petrovic et al. [91], Hernández et al.
23 [92], Farré et al, [93] and Xu et al. [51].

1 3.2.1 Ionization techniques in mass spectrometry

2 There are several techniques by which ions can be generated and mass analysis can be performed. The most
3 commonly used technique for the identification of POPs is electron ionization (EI), which is typically
4 coupled with GC. EI mass spectral libraries have been developed for matching spectra for non-target
5 compound identification. Significant fragmentation occurs under EI conditions, which is crucial to the
6 success of library matching. However, it can also result in elimination of the molecular ion and thus, the
7 possibility of obtaining an elemental composition. Increasing the pressure in the ion source leads to
8 chemical ionization (CI), which is a much softer form of ionization that can result in an improved yield of
9 molecular ions as well as, in some cases, diagnostic ion-molecule reactions [94]. Electron capture negative
10 chemical ionization (ENCI) is a particularly helpful technique for the determination of halogenated POPs
11 [94, 95] because it is selective towards compounds with a high electron affinity. Despite the selectivity
12 benefit, ENCI is not usually employed for routine methods and dissociative electron capture into Cl⁻ and
13 Br⁻ can hamper analysis. Other forms of “softer” ionisation include field ionisation (FI) and photo ionisation
14 (PI). Field ionisation uses high electrical fields in the close proximity of emitter needles. During
15 photoionisation molecules are ionised by absorption of two or more ultraviolet (UV) or one single vacuum
16 ultraviolet (VUV) photon. These forms of ionisation make detectors highly sensitive and selective for
17 aromatic species, however the laser based photo-ionisation technologies are still too expensive and
18 sophisticated for a broader application [96].

19 Liquid chromatography (LC) and the associated atmospheric ionization methods such as electrospray (ESI)
20 are now widely used for the determination of polar non-volatile and/or thermally labile POPs. Additional
21 ionisation techniques include atmospheric pressure chemical ionisation (APCI), and atmospheric pressure
22 photo ionisation (APPI). APPI is typically used for non-polar compounds like halogenated flame retardants
23 and PAH. APCI and APPI have mainly been coupled to LC, but recently these ionization techniques are
24 also being used for GC analysis [97]. Methods have been developed for specific POPs, including dioxins
25 [98, 99] and polybrominated diphenyl ethers [100, 101].

3.2.2 Nominal mass resolution mass spectrometers

Mass analysers are broadly classified as being capable of performing either nominal or high mass resolution measurements. The resolution of a mass spectrometer ($R=M/\Delta M$) is defined by the degree of peak overlap at a 10% valley (10%V) or full width of a peak at half peak height (full width half maximum - FWHM). Quadrupole and quadrupole ion trap mass spectrometers usually provide nominal mass resolution ($R = 100 - 1000$). High resolution mass spectrometers include magnetic sector (HRMS) and some time-of-flight (HRTOF) instruments, which are capable of operating at $R=10,000 - 50,000$. Fourier transform ion cyclotron resonance (FTICR) and Orbitrap mass spectrometers may be considered *ultrahigh* resolution instruments ($R = 100,000 - 1,000,000$), which can resolve mass spectral peaks that differ by as little as 1 mDa [102].

Mass resolution is important but it is not the only parameter that needs to be considered. The sensitivity of the instrument is also important as there is no benefit to being able to distinguish two compounds of a similar mass, if you are unable to detect them in a sample at environmentally relevant concentrations. The linear range of a detector is also important, in environmental investigations there are often several orders of magnitude difference in concentrations between samples and also between individual congeners of the same compounds class. Whilst some of these issues can be overcome by dilution or concentration and reanalysis, there are obvious time and cost benefits of not having to incur this extra level of sample preparation and analysis.

Quadrupole instruments were the first bench-top instruments, becoming commercially available in the 1980s, which significantly reduced the cost of mass spectrometry detection. They are now very wide spread and can be found in most commercial and research laboratories. A quadrupole instrument is typically limited to unit mass resolution measurements, but its excellent linear dynamic range (4-6 orders of magnitude), sub picogram sensitivity, ease of use and low operating costs make it a popular choice for initial screening using targeted compound analysis and full scan applications. They can also be a cost

1 effective option when analysing samples where low detection limits are not essential, such as the
2 determination of PBDEs [103] and PAHs [104] in dust.

3 Tandem instruments usually consist of two quadrupole mass analysers, sandwiching an ion guide collision
4 cell. Other hybrid types incorporate a quadrupole ion trap, time-of-flight, FT-ICR or Orbitrap as the final
5 mass analyser. The first quadrupole is used for initial selection and focusing. The collision cell is filled with
6 a collision gas (helium, nitrogen, oxygen or argon) to produce specific fragments that are subsequently
7 mass separated by the final mass analyser. Tandem mass spectrometry has the advantage of being able to
8 filter out non target ions and can therefore reach low femtogram detection levels because the chemical noise
9 is substantially reduced. Most MS/MS instruments exhibit a good linear dynamic range (4-6 orders of
10 magnitude) which makes this technique an excellent choice for targeted determination of trace level
11 contaminants; such as PBDEs and other POPs in biota and food webs [105-107]. Drawbacks of MS/MS
12 instruments are that MS/MS experiments are more complex to develop and optimize than single stage MS
13 experiments, and non-target compounds are not detected.

14 3.2.3 High resolution mass spectrometers

15 The magnetic deflection instruments are capable of both HRMS and tandem MS, but usually only the
16 former mode is used for targeted compound analysis in environmental and environmental forensics
17 applications. Modern magnetic sector instruments are routinely tuned to achieve a resolution of >10,000
18 (10% valley), and they can achieve low-femtogram instrument detection limits (IDLs) or even attogram
19 IDLs when coupled to cryogenic zone compression gas chromatography [108]. Such instruments also
20 exhibit an excellent dynamic range of 6 to 7 orders of magnitude. However, most analysis is undertaken
21 using selective ion monitoring (SIM) where the detector only focuses on several target ions. This has helped
22 to establish HRMS as the technique specified by the U.S. EPA for low level congener specific analysis of
23 compounds such as PCBs (Method 1668C), PCDD/F (Method 1613) and PBDEs (Method 1614). Whilst
24 SIM increases the sensitivity by a factor of 10 to 100 it does mean that non-target compounds are not

1 detected. Nowadays, the main drawbacks of magnetic deflection instruments are commercial availability
2 and the significant costs associated with maintenance and user training.

3 TOF instruments are relatively simple in design and easy to maintain, tune, calibrate and operate. They are
4 available as low resolution, mid resolution (7,000 to 10,000) and high resolution ($> 20,000$) instruments.
5 Recent improvements over the last decade have allowed low femtogram detection levels to be achieved.
6 The major disadvantage of most TOF instruments is the low dynamic range which is usually limited to 4
7 orders of magnitude. The chief advantage of the TOF is its high data acquisition rate (20-200 scans/second)
8 while recording full-scan mass spectra. This makes them the detector of choice for fast GC and
9 environmental forensics GCxGC applications for POPs [31]. The ability of the TOF to collect data for an
10 entire mass range gives them the excellent potential as a tool for non-targeted screening and for the
11 simultaneous determination of multiple POPs [82, 107, 109, 110].

12 There are two instrument types available that can operate at ultra-high resolution. These include Fourier
13 transform ion cyclotron resonance (FT-ICR) and orbitrap mass analyzers. Both are capable of achieving a
14 mass resolution of $>100,000$. Newer instruments are quite sensitive, with LC and GC amenable compounds
15 detectable in the mid femtogram range. Since these detectors work by trapping a finite number of ions in a
16 cell, they suffer from space charge effects which limit the instrument to approximately 4 orders of
17 magnitude. The acquisition rate and resolving power of these instruments are limited by the length of the
18 transient, i.e., the time the ions are confined during detection. Longer transients result in higher resolution.
19 The timescale is most compatible with LC and direct infusion (> 1 scan/second); but GC experiments are
20 still possible. For example Peterson et al. [111] used a GC-Enabled QLT-Orbitrap for the determination of
21 dioxins in environmental samples. Due to their exceptional high mass resolution FT-ICR and orbitrap mass
22 analyzers instruments are well suited to the determination of complex samples or when an unknown sample
23 needs to be characterised [107, 109, 110].

1 3.2.4 Comparison of different detectors for the determination of POPs

2 There is not one technique that can be considered the best for the forensic determination of POPs, the
3 various advantages of the different techniques are presented above and have been summarised in Table 4.
4 Whilst the extraction method used will have a large influence on the method detection limits a summary of
5 recent detection limits in various environmental matrices for different POPs is included. Reviews on the
6 determination of POPs have been undertaken by many authors, recent examples include; Reiner et al. [9],
7 Xu et al. [51], Farre et al. [93], Tang [50]. The vast majority of methods reported in the analytical literature
8 have been developed to meet a set of criteria constant with a specific regulatory limit. The method
9 development is therefore focused on proving that the technique has detection levels below this limit so there
10 can be confidence if a sample has exceeded it or not. Whilst this makes sense for a screening exercise it
11 does mean that some diagnostic information which would be useful for a forensic investigation may not be
12 recorded as much lower detection limits could be achieved using the same instrument. A good example of
13 this is for the determination of pesticide residues in foods. The majority of methods on different instruments
14 all present a limit of detection in the range of <10 ng/g [57, 112-117] to coincide with regulatory maximum
15 residue limits. However, there are examples where much lower detection limits of <500 pg/g can be
16 achieved for these compounds [86, 118].

17 As well as focusing on a specific regulatory limit, many methods focus on one compound class. This type
18 of measurement favours techniques such as MS/MS which can filter the data to get excellent signal to noise
19 ratios and low detection limits for a specific set of compounds. A recent example of this was Organtini et
20 al. [99] who used APGC-MS/MS to detect 2 fg of dioxins (on column) resulting in detection limits of <1
21 ng/kg in soil. Again whilst this makes sense for a screening exercise it does mean that some diagnostic
22 information which would be useful for a forensic investigation may not be recorded. There are several
23 advantages of collecting full scan data on an instrument such as TOFMS. This can allow for the
24 determination of multiple POPs in the same run [80]. Data can also be re-evaluated later to check for the
25 presence of other potential contaminants of concern without the need for re-analysis. This allowed Megson

1 et al. [119] to identify an additional group of toxic organophosphates in jet oil following an initial screening
2 exercise in a forensic investigation. Finally if a new potential POP is discovered its presence can be checked
3 for in historic samples to identify when it arrived in the environment. Nowadays some instruments are able
4 to combine multiple functions and collect both MS/MS data and full scan TOF data, this may prove to be a
5 very powerful tool in future investigations.

6

1 **Table 4. Comparison of different detectors used for the determination of POPs**

		Detector type					
		ECD	qMS	TOFMS	HRMS (magnetic deflection)	MS/MS	Orbitrap
Mass resolution		na	nominal	nominal → >20,000	>10,000	nominal → >20,000	> 100,000
Linier range (orders of magnitude)		4 to 5	4 to 6	3 to 4	6 to 7	4 to 6	3 to 4
Cost of analysis		low	low	low → high	high	high	high
Selectivity		v.low	low	high	high	high	high
Sensitivity		high (halogenated) → low (non-halogenated)	moderate	moderate → high	high	high	moderate → high (depending on acquisition settings)
Limit of detection	OC Pesticides, including; Chlordane, Dieldrin, Endrin, Heptachlor, Hexachlorobenzene, Mirex, Toxaphene, DDT, α-Hexachlorocyclohexane, β-Hexachlorocyclohexane, γ- Hexachlorocyclohexane (Lindane), Chlordecone, Pentachlorobenzene, Endosulfan	<10 ng/g (various foods) [57,112,113] <5 pg/g [118], <500 pg/g* [86] (soil) 8-33 pg/g (fish) [120]	<10 ng/g (various foods) [114, 115] <500 pg/g (soil) [118] <1 ng/g (soil) [121]	<50 pg / L (water) [122]	<10 ng/L (serum) [45, 123]	<10 ng/g (various foods) [113, 115, 116]	<2 µg/L (analytical standard) [124] 10 ng/g (various foods) [117]
	Polychlorinated biphenyls (PCBs) and **Polychlorinated naphthalenes (PCNs)	<140 pg/g, [125] (sediment) <1.6 ng/g* [86] (soil) <2.4 µg/L (water) [126] 10 pg/g (fish) [127] 8-33 pg/g (fish) [120]	<1 ng/g (soil) [121] < 0.14 µg/L (water) [126]	1 to 50 pg/µL * (analytical standard) [32] 0.5 to 10 pg/µL *** (analytical standard) [82] 1-15 pg/µL *** (serum) [82]	<80 ng/g (soil) [128] <30 pg/L (water) [128] <50 ng/L (serum) [45, 123]	< 3 pg/g (sediment) [125]	only reported LOD of 10 fg for octafluoronaphthalene [111]

Polychlorinated dibenzo-p-dioxins and **polychlorinated dibenzofurans,	< 5.2 pg/g* (sediment) [89]	< 1 pg/g fat (beef) [129]	<0.5 pg/μL * (2,3,7,8-TCDD analytical standard) [130]	0.54-5.04 ng/g (sediment)[131] <2 ng/L (serum) [123]	2 fg/ul (analytical standard) & <1 pg/g (soil) [99]	
Polybrominated diphenyl ethers (PBDEs) and Hexabromobiphenyl	8-33 pg/g (fish) [120]	0.005-0.1 μg/L (water) [132] 0.25-5 ng/g (soil) [132] 10 pg/g (fish) [127]	0.025-5 ng/g * (fish) [133]	<200 ng/g (soil) [128] <50 pg/L (water) [128] <10 ng/L (serum) [123]	<0.1 ng/g (sediment) [134] 1-80 ng/g (dust) [135] 1-25 ng/g (dust) [136]	1-250 pg/g (fish) [137]
Hexabromocyclododecane (HBCD)	6 pg/g (milk) [138]	0.03 μg/L (water) [139]	<30 pg/g**** (fish) [61] 0.9-4.5 pg/g (fish) [61]	5 ng/g (dust) [139]	<9 pg/g (fish) [61] 0.9-2.1 pg/g (fish) [61] 1 ng/g (dust) [135] 1-25 ng/g (dust) [136]	10 pg/g (fish) [137] <13 pg/g (fish) [61] 1.1-3.0 pg/g (fish) [61]
**Perfluorooctanesulfonic acid (PFOS), and its salts **perfluorooctanesulfonyl fluoride (PFOSF)	Response detected – no LOD reported [51]	0.05 μg/L (water) [141]	0.027-2.820 ng/g (dust) [142]		0.01–2.00 ng/L (water) [143] 0.04–8.00 ng/g (sediment) [143] 0.02–2.26 ng/g (biota) [143] 5-50 pg/g (fish) [67] 0.25-1 ng/g (dust) [136]	0.3-3 ng/L (water) [144] 5 to 200 ng/L [68]

- 1 *analysis by GCxGC-μECD
- 2 ** reported values for one congener, not total of all congeners
- 3 ***analysis by GCxGC-TOFMS
- 4 ****analysis by HRTOFMS
- 5

1 4 Discussion

2 4.1 Method development, validation and QA/QC

3 In environmental forensics investigations it is especially important to follow strict QA/QC procedures and
4 keep records of the work undertaken. In forensics investigations, a laboratory is required to have a
5 comprehensive quality assurance / quality control (QA/QC) program that often includes accreditation to
6 ISO 17025. Accreditation ensures that laboratories are following proper procedures and participating in
7 performance evaluation studies to assess data quality. Whilst this makes perfect sense from a legal
8 standpoint and in routine screening and compliance, it may not always be achievable in a forensics
9 investigation. There may be instances where new and unknown chemicals have been detected, in these
10 cases there may not be analytical standards or an accredited method available for analysis.

11 Methods are characterized by their data quality objects which include sensitivity (detection limits), and
12 selectivity (precision or uncertainty and accuracy). Analytical cost and speed of analysis are also important
13 considerations when selecting or developing a method. Just because a method is accredited it does not
14 specifically mean it is the method that should be selected. For example if a lab is accredited to undertake
15 EPA Method 608 (GC-ECD) or 625 (GC-MS) for PCBs this may be suitable for a routine screening.
16 However, in a forensics investigation the methods are unlikely to be able to provide the detection limits
17 required, may underestimate actual PCB concentrations, would not necessarily capture degraded PCBs that
18 contain non-Aroclor congeners, and may therefore lead to erroneous quantification of Aroclors that are not
19 necessarily attributed to the source of contamination [145].

20 In some cases analytical methods can be overlooked because they haven't been accredited. This can be
21 detrimental if the method is capable of providing better quality, more useful data than the accredited
22 method. However, if a non-accredited method is used can the court still have the same confidence in the
23 quality of the data or applicability of the method? In 1993, the United States Supreme Court issued a
24 landmark decision in *Daubert v. Merrell Dow Pharmaceuticals, Inc.*, 509 United States (US) 579, which

1 made significant changes in the standards for admissibility of expert opinions in federal courts. These
2 changes included a gatekeeping requirement in Rule 104(a) under which courts must screen expert opinions
3 for reliability and exclude ‘junk science.’ In Daubert, the Supreme Court also established a new, more
4 flexible set of criteria for reliability and admissibility of expert opinion [146-149].

5 According to Kanner [148], for a method to be accepted in court it needs to satisfy the following Daubert
6 criteria:

- 7 • the scientific technique (or “theory” applied) must be testable and verifiable;
- 8 • the technique has been published in a peer-viewed journal or other similar publication;
- 9 • the technique has a defined rate or margin of error;
- 10 • the technique is used with appropriate standards and controls;
- 11 • and, the scientific community has accepted the technique or theory to a significant degree.

12 In all cases the production of accurate and precise results is critical and a laboratory must have a
13 comprehensive quality assurance / quality control (QC/QC) program. Accreditation of a method is
14 preferable, however non-accredited methods should not be discounted as they have and will continue to be
15 accepted in a court of law. The most important requirement for any method is that it is fit for the purpose it
16 was intended i.e. it satisfies all the data quality objectives needed for the client to make correct decisions
17 using the data produced by the method [3]. Different jurisdictions allow different types of methods,
18 procedures and instrumentation for regulatory work. Some jurisdictions (e.g. United States) typically
19 require prescriptive methods and others (e.g. Ontario, Canada) allow more flexible methods and procedures
20 that are performance based, as long as they are fully validated. Method validation requires a number of
21 steps that involve the testing of spiked matrix samples and certified reference materials to determine method
22 characteristics (detection limits, accuracy, and uncertainty) and to ensure the method is rugged, generally
23 free of interferences and data are reproducible. Details on the various steps for method development and
24 validation can be obtained in Reiner et al [9] and HimaBindu et al [150].

1 4.2 PCBs as an example for environmental forensics investigations

2 PCBs are one of the most widely studied group of POPs; however there is no commonly used method for
3 PCB determination. PCBs were first detected as interferences in organochlorine pesticide GC-ECD
4 chromatograms and it took almost 7 years to identify and confirm PCBs as environmental contaminants
5 [9]. There are such a large number of PCB methods currently in use that when comparing data it is rare to
6 find laboratories that used the same method. Early methods such as EPA Method 608 or 8082A used GC-
7 ECD with PCB technical mixtures like Aroclors to calibrate and quantify results [44], however there are
8 many limitations to these methods, which are discussed in Johnson et al. [145]. There are relatively few
9 methods, such as EPA Method 1668C, that use GC-MS and GC-HRMS and are designed to determine all
10 209 PCBs, instead the vast majority of congener methods determine a reduced set of congeners. One of the
11 main justifications for this has been because all 209 PCBs are not found in detectable concentrations in
12 environmental samples. Out of the 5 Aroclors (A1016, A1242, A1248, A1254 and A1260) that contribute
13 to >97% of Monsanto's sales [151] only 157 of the 209 possible PCBs are present, and only about 130-150
14 PCBs have ever been detected in environmental samples [9, 31].

15 Many commercial laboratories provide limits of detection for PCBs in the range of 10 – 1000 $\mu\text{g kg}^{-1}$. This
16 can be a significant limitation as background total PCB concentrations (in UK urban soils) are between 0.01
17 – 40 $\mu\text{g kg}^{-1}$ [152]. Therefore, sample clean-up and analysis by HRMS is often required to improve limits
18 of detection to less than 1 $\mu\text{g kg}^{-1}$ by removing or filtering out many interfering compounds. In many
19 environmental forensics investigations there can be multiple sources of contamination which makes source
20 identification a challenge, especially if the original signature has been altered by post depositional processes
21 such as; volatilization and dispersion, biodegradation, uptake in biota, and biotransformation and
22 elimination. In such cases it is important to select a method that has the required sensitivity but is also able
23 to quantify as many individual congeners as possible to improve the potential power of the statistical
24 methods used. In this instance GCxGC is a very useful tool as it provides an extra dimension of separation
25 which significantly increases the resolving capacity. Using GCxGC has allowed for the identification of

1 over 190 of the 209 PCB congeners along with simultaneous identification of other organohalogenated
2 contaminants [81, 31]. This method has been successfully applied to several environmental forensics
3 applications such as determining the provenance of sea birds [87] and establishing the source of PCB
4 exposure in humans [32].

5 4.3 Determination of Polychlorinated Dioxins and Furans

6 Although dioxins, furans and PCBs are chemically similar, most dioxin and furan methods are generally
7 much more accurate and precise than most PCB methods. The majority of dioxin and furan methods are
8 based on US EPA 1613 and use isotope dilution with GC-HRMS to determine the seventeen 2,3,7,8-
9 substituted toxic congeners. Because of the highly toxic nature of these compounds, they are regulated
10 globally and are often detected in exposure and contamination assessments and environmental forensics
11 investigations [152-156].

12 Due to the extensive sample concentration and cleanup steps required for dioxin methods, method
13 turnaround times are long and costs are significantly higher than those for other types of analytes [41, 42].
14 Detailed contaminated site characterization for reclamation can require a significantly large number of
15 samples to delineate the contaminated area, which can be extremely costly. Bioassays have been used, but
16 can require extensive sample extract cleanup to produce accurate results. The QuEChERS procedure uses
17 a modified extraction procedure to significantly reduce sample preparation time. Haimovici et al. [131];
18 developed a method for site cleanup that significantly enhances productivity and reduces analytical
19 turnaround time by a factor of 3 or more [154, 131].


20 4.4 Non target-analysis and identification of POPs

21 There are a few hundred chemicals that make up the Stockholm POPs Convention list (Table 1) [1]. These
22 compounds are typically determined using targeted methods [90]. Using targeted analysis, many thousands
23 of compounds are overlooked; this may be detrimental in an environmental forensics investigation where
24 the source of contamination is unknown. The development of time of flight mass spectrometers [109] and

1 trapping instruments (Orbitrap) [157] to increase sensitivity, scan speed and resolution enables them to be
2 used for non-target full scan analysis. The use of scanning instruments (TOF, Q-TOF, Orbitrap) enables
3 analysts to obtain full scan spectra which can be reviewed manually to obtain unique mass to charge ratios
4 and assign elemental compositions for the suspect compounds being investigated [109, 157-160].
5 Comparing the elemental compositions of the fragment ions to the molecular ion and searching databases
6 such as Chemspider (<http://www.chemspider.com/>) enables the postulation of a possible structure. This can
7 be used for the discovery of new contaminants of concern and may also identify compounds that may be
8 unique to a specific source. There have been many instances where chemical impurities or co-contaminants
9 present in samples have revealed some key information in an environmental forensics investigation. Recent
10 examples include 2,4,6,8-Tetrachlorodibenzothiophene as an impurity in 2,4-D production from a specific
11 manufacturer in Newark Bay [161] impurities in dechlorane plus production from a manufacturer in the
12 Niagara region [162], and xylenyl cresyl phosphates in used jet oil [119]. Had just the main contaminant in
13 question been determined in these investigations then this key diagnostic information would not have been
14 recorded. There are also an increasing number of studies that have used non-targeted methods to identify
15 mixed brominated / chlorinated aromatic compounds such as dioxins, furans and PAHs in environmental
16 samples [107, 110]. Using complimentary non-targeted and targeted methods can help to establish more
17 accurate levels of risk and provide more diagnostic information for forensic investigations.

18 Mass spectra generated from complex samples can contain many thousand fragment ions. Mass defect plots
19 (MDPs) can be used to process this data and identified many compounds in a single analysis. The mass
20 defect is the difference between the nominal and exact mass of a chemical element or compound. Plotting
21 the mass defect against the mass to charge ratio allows similar compounds (homologues and congeners) to
22 be identified [163]. MDPs can be calibrated to a variety of elements (CH_2 , Cl, Br, CF_2) as shown in Figure
23 2 [107, 164-169]. Once one of the compounds in a row or column is identified, the remaining compounds
24 of each connected row or column can relatively easily be identified as they are likely to be homologues.

1 These techniques are very useful when samples are chemically complex or when the main compounds of
2 interest have not yet been established.

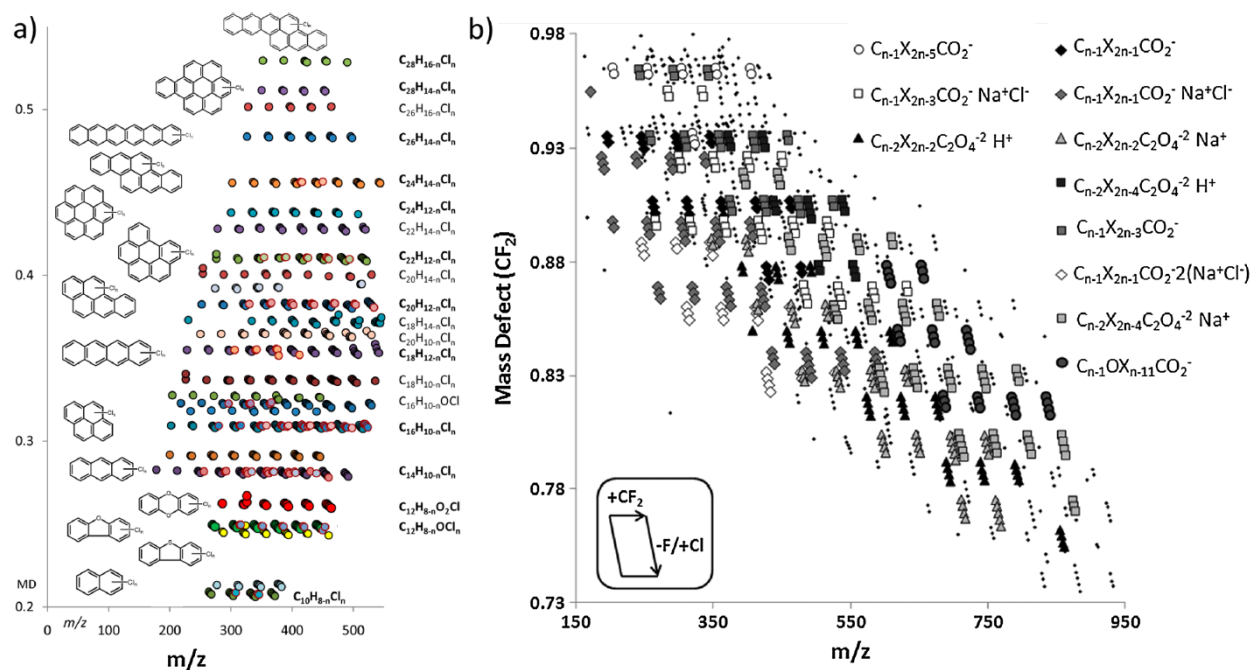

3
4 **Figure 2.** Mass Defect plot using the H/Cl scale of an extract from a plastic manufacturing plant fire (a),
5 and using the CF₂ scale in a perfluorinated compounds thermolysis study (b) Reproduced with permission
6 from Ortiz et al. [109] (original figures a) from Fernando et al. [110] and b) from Myers et al. [169])

7 5 Conclusions

8 The determination of persistent organic pollutants for environmental forensics investigations is not
9 straightforward. There are a wide variety of accredited and non-accredited methods available and it can be
10 difficult deciding which one to choose. There is not one analytical technique or method that can solve all
11 the questions asked in an investigation. The decision needs to be made on a case specific basis depending
12 on the required detection limits, linier range, sensitivity and selectivity. Unfortunately there are times when
13 this decision is made by those with little understanding of the different techniques and their decision is
14 based primarily on the cost of the analysis. Taking the time to plan an investigation properly and understand

1 the merits of the different techniques, will lead to the generation of better quality data which will be more
2 likely to uncover the truth and may also turn out to be the more cost effective option.

3

4 References

5 [1] Stockholm Convention Secretariat, (2001), UNEP, <http://chm.pops.int/> , Accessed Sept 2015

6 [2] Howard, PH. Muir, DCG. (2010) Identifying New Persistent and Bioaccumulative Organics Among
7 Chemicals in Commerce. *Environ. Sci. Technol.* 44 (7):2277-2285

8 [3] Bethem, R. Boison, J. Gale, J. Heller, D. Lehotay, S. Loo, J. Musser, S. Price, P. Stein, S. (2003)
9 Establishing the fitness for purpose of mass spectrometric methods. *J. Am. Soc. Mass Spectrom.*
10 14: 528-541

11 [4] Mudge, SM. (2008) *Methods in Environmental Forensics*. CRC Press

12 [5] Morrison, RD. Murphy, BL. (2006) *Environmental Forensics: Contaminant Specific Guide*. Academic
13 Press

14 [6] Murphy, BL. Morrison, RD. (2014) *Introduction to Environmental Forensics*, 3rd Edition. Elsevier

15 [7] Simmons, B. (1999) Are the data legally defensible? *Hydrovisions*. 8 (2), 8-9

16 [8] Wait, D. Ramsey, C. Maney, J. (2015) The Measurement Process. In *Introduction to Environmental*
17 *Forensics (3rd Edition)* Morrison, RD. Murphy, BL. (eds). Elsevier.

18 [9] Reiner, EJ. Jobst, KL. Megson, D. Dorman, FL. Focant, JF. (2013) Chapter 3, *Environmental Forensics*
19 *Analytical Methodology of POPs*. In: *Environmental Forensics for Persistent Organic Pollutants*.
20 Elsevier

21 [10] Poole, CF. (2003) New trends in solid-phase extraction, *TrAC, Trends. Anal. Chem* 22: 362-373

22 [11] Pujadas, E. Diaz-Ferrero, J. Marti, R. Broto-puig, F. Comellas, L. Rodriguez-Larena, MC. (2001)
23 Application of the new C18 speedisks™ to the analysis of polychlorinated dibenzo-p-dioxins and
24 dibenzofurans in water and effluent sample. *Chemosphere*. 43: 449-454

- 1 [12] Focant, JF. DePauw, E. (2002) Fast automated extraction and clean-up of biological fluids for
2 polychlorinated dibenzo-p-dioxins, dibenzofurans and coplanar polychlorinated biphenyls
3 analysis. *J. Chromatogr. B.* 776: 199-212
- 4 [13] Prieto, A. Zuloaga, O. Usobiaga, A. Etxebarria, N. Fernandez, LA. (2007) Development of a stir bar
5 sorptive extraction and thermal desorption-gas-chromatography-mass spectrometry method for the
6 simultaneous determination of several persistent organic pollutants in water samples. *J.*
7 *Chromatogr. A.* 1174: 49-49.
- 8 [14] Li, G. Zhang, L. Zhang, Z. (2008) Determination of polychlorinated biphenyls in water using dynamic
9 hollow fiber liquid-phase microextraction and gas chromatography-mass spectrometry.
10 *J.Chromatogr. A.* 1204: 119-122
- 11 [15] Arthur, CL. Pawliszyn, J. (1990) Solid phase microextraction with thermal desorption using fused
12 silica optical fibres, *Anal. Chem.* 62: 2145-2148
- 13 [16] Lohmann, R. (2012) A critical review of low-density polyethylene's partitioning and diffusion
14 coefficients for trace organic contaminants and implications for its use as a passive sampler.
15 *Environ. Sci. Technol.* 46, 606-618.
- 16 [17] Muir, D. and Sverko, E. (2006). Analytical methods for PCBs and organochlorine pesticides in
17 environmental monitoring and surveillance: a critical appraisal. *Anal Bioanal Chem* 386:769–789
- 18 [18] Melymuk, L. Robson, M.. Helm, P.A. Diamond, M.L. (2012) PCBs, PBDEs, and PAHs in Toronto
19 air: spatial and seasonal trends and implications for contaminant transport. *Science of the Total*
20 *Environment* 429: 272-280.
- 21 [19] Khairy, M.A. and Lohmann, R. (2013) Feasibility of using low density polyethylene sheets to detect
22 atmospheric organochlorine pesticides in Alexandria, Egypt. *Environmental Pollution.* 181: 151-
23 158.
- 24 [20] Harner, T. Shoeib, M. Diamond, M.L. Ikonou, M. Stern, G. (2006) Passive sampler derived air
25 concentrations of PBDEs along an urban–rural transect: spatial and temporal trends. *Chemosphere.*
26 64(2): 262-267.

- 1 [21] Meire, R.O. Khairy, M. Targino, A.C. Galvao, P.M.A. Torres, J.P.M. Malm, O. Lohmann, R (2016).
2 Use of passive samplers to detect organochlorine pesticides in air and water at wetland mountain
3 region sites (S-SE Brazil) *Chemosphere* 144. 2175-2182
- 4 [22] Khairy, M. Muir, D. Teixeira, C. Lohmann, R. (2014). Spatial Trends, Sources, and Air-Water
5 Exchange of Organochlorine Pesticides in the Great Lakes Basin Using Low Density Polyethylene
6 Passive Samplers. *Environmental Science & Technology*. 48(16). 9315-9324
- 7 [23] Antunes, P. Viana, P. Vinhas, T. Capelo, JL. Rivera, J. Gaspar, EMSM, (2008) Optimization of
8 pressurized liquid extraction (PLE) of dioxins-furans and dioxin-like PCBs from environmental
9 samples. *Talanta*. 75: 916-925
- 10 [24] Saito, K. Sjodin, A. Sandau, CD. Davis, MD. Nakazawa, H. Matsuki, Y. Patterson Jr, DG. (2004)
11 Development of a accelerated solvent extraction and gel permeation chromatography analytical
12 method for measuring persistent halogenated compounds in adipose and organ tissues.
13 *Chemosphere*. 57: 373-381
- 14 [25] Richter, BE. Jones, BA. Ezzell, JL. Porter, NL. Avdalovic, N. Pohl, C. (1996) Accelerated solvent
15 extraction: A technique for sample preparation. *Anal. Chem*. 68: 1033-1039
- 16 [26] Srogi, K. (2007) A Review: application of microwave techniques for environmental analytical
17 chemistry, *Anal. Lett*. 39: 1261-1288
- 18 [27] Sun, Y. Takaoka, M. Takeda, N. Matsumoto, T. Oshita, K. (2005) Application of water as a solvent in
19 microwave-assisted extraction for analysis of PCBs and CBzs in fly ash. *J. Sep Sci*. 28: 585-588
- 20 [28] Vazquez-Roig, Pablo, and Yolanda Picó. "Pressurized liquid extraction of organic contaminants in
21 environmental and food samples." *TrAC Trends in Analytical Chemistry* 71 (2015): 55-64.
- 22 [29] Miyawaki, T. Kawashima, A. Honda, K. (2008) Development of supercritical carbon dioxide
23 extraction with a solid phase trap for dioxins in soils and sediments. *Chemosphere*. 70: 648-655.
- 24 [30] Garcia-Rodriguez, D. Carro-Diaz, AM. Lorenzo-Ferreira, RA (2008) Supercritical fluid extraction of
25 polyhalogenated pollutants from aquaculture and marine environmental samples: A review. *J. Sep.*
26 *Sci*. 31: 1333-1345

- 1 [31] Megson, D. Kalin, RB. Worsfold, PJ. Gauchotte-lindsay, C. Patterson Jr, DG. Lohan, MC. Comber, S.
2 Brown, TA. O’Sullivan, G. (2013) Fingerprinting polychlorinated biphenyls in environmental
3 samples using comprehensive two-dimensional gas chromatography with time-of-flight mass
4 spectrometry. *J. Chromatogr. A.* 1318: 276-283.
- 5 [32] Megson, D. Focant, J.F. Patterson Jr, DG. Robson, M. Lohan, MC. Worsfold, PJ. Comber, S. Kalin,
6 R. Reiner, E. O’Sullivan, G. (2015) Can polychlorinated biphenyl (PCB) signatures and enantiomer
7 fractions be used for source identification and to age date occupational exposure? *Environ. Int.* 81:
8 56–63
- 9 [33] Anastassiades M, Lehotay SJ, Stajnbaher D, Schenck FJ. 2003 Fast and easy multiresidue method
10 employing acetonitrile extraction/ partitioning and “dispersive solid-phase extraction” for the
11 determination of pesticide residues in produce. *JAOAC Int* 86:412–431
- 12 [34] Pinto, CG. Laespada, MEF. Martín, SH. Ferreira, AMC. Pavón, JLP. Cordero, BM. (2010) Simplified
13 QuEChERS approach for the extraction of chlorinated compounds from soil samples. *Talanta.* 81:
14 385-391
- 15 [35] Norli, HR. Christiansen, A. Deribe, E. (2011) Application of QuEChERS method for extraction of
16 selected persistent organic pollutants in fish tissue and analysis by gas chromatography mass
17 spectrometry. *J. Chromatogr. A.* 1218: 7234-7241
- 18 [36] Correia-Sá, L. Fernandes, VC. Carvalho, M. Calhau, C. Domingues, VMF., Delerue-Matos, C. (2012)
19 Optimization of QuEChERS method for the analysis of organochlorine pesticides in soils with
20 diverse organic matter. *J. Sep. Sci.* 35: 1521-1530
- 21 [37] Bruzzoniti, MC., Checchini, L., De Carlo, RM., Orlandini, S., Rivoira, L., Del Bubba, M., (2014).
22 QuEChERS sample preparation for the determination of pesticides and other organic residues in
23 environmental matrices: a critical review. *Analytical and Bioanalytical Chemistry.* 406, 4089-4116
- 24 [38] Gonzalez-Curbelo, MA., Socas-Rodriguez, B., Herrera-Herrera, AV., Gonzalez-Salmo, J., Hernandez-
25 Borges, J., Rodriguez-Delgado, MA., 2015. Evolution and applications of the QuEChERS method.
26 *Trends in Analytical Chemistry.* 71, 169-185.

- 1 [39] Poole, SK. Dean, TA. Oudsema, JW. Poole, CF. (1990) Sample preparation for chromatographic
2 separations: an overview. *Anal. Chim. Acta* 236: 3-42
- 3 [40] Van Leeuwen, SPJ. de Boer, J. (2008) Advances in the gas chromatographic determination of persistent
4 organic pollutants in the aquatic environment. *J. Chromatogr. A.* 1186: 116-182.
- 5 [41] Reiner, EJ. Clement, RE. Okey, AB. Marvin, CH. (2006) Advances in the analytical techniques for
6 polychlorinated dibenzo-p-dioxins, polychlorinated dibenzofurans and dioxin-like PCBs. *Anal.*
7 *Bioanal. Chem.* 386: 791-806.
- 8 [42] Reiner, EJ. (2010) Analysis of dioxin and related compounds, *Mass Spectrometry Reviews.* 29: 526-
9 559
- 10 [43] Muir, D. Sverko, E. (2006) Analytical methods for PCBs and organochlorine pesticides in
11 environmental monitoring and surveillance: A critical appraisal. *Anal. Bioanal. Chem.* 386: 769-
12 789.
- 13 [44] Clement, RE. Reiner, EJ. Bhavsar, SP. (2012) Organohalogen contaminants of emerging concern in
14 Great Lakes fish: A review. *Anal. Bioanal. Chem.* 404: 2639-2658
- 15 [45] Wittsiepe, J., Nestola, M., Kohne, M., Zinn, P., Wilhelm, M. (2014). Determination of polychlorinated
16 biphenyls and organochlorine pesticides in small volumes of human blood by high-throughput on-
17 line SPE-LVI-GC-HRMS. *Journal of Chromatography B.*
- 18 [46] Kole, PL., Venkatesh, G., Kotecha, J., Sheshala, R. (2011). Recent advances in sample preparation
19 techniques for effective bioanalytical methods. *Biomedical Chromatography*, 25(1-2), pp.199-217
- 20 [47] Lambropoulou, DA. (2010) An overview of modern extraction techniques for the determination of
21 organic pollutants in environmental matrices: a Review, *Current Organic Chemistry.* 14: 2247-
22 2267
- 23 [48] Fidalgo-Used, N. Blanco-Gonzalez, E. Sanz-Medel, A. (2007) Sample handling strategies for the
24 determination of persistent organic pollutants from biota samples. *Anal. Chim. Acta.* 590: 1-16
- 25 [49] Rawa-Adkonis, M. Wolska, L. Namiesnik, J. (2003) Modern techniques of extraction of organic
26 analytes from environmental matrices. *Crit. Rev. Anal. Chem.* 33(3): 199-248

- 1 [50] Tang, HPO., (2013). Recent development in analysis of persistent organic pollutants under the
2 Stockholm Convention. *Trends in Analytical Chemistry*, 45, 48-66.
- 3 [51] Xu, W., Wang, X., Cai, Z., (2013). Analytical chemistry of the persistent organic pollutants identified
4 in the Stockholm Convention: A review. *Analytica Chimica Acta*. 790, 1-13
- 5 [52] Zwiener, C. Frimmel, FH. (2004) LC-MS analysis in the aquatic environment and in water treatment
6 – a critical review. Part I: Instrumentation and general aspects of analysis and detection. *Anal.*
7 *Bioanal. Chem.* 378; 851-861
- 8 [53] Jaume, A. Stampachiachiere, S. Pérez, S. Barceló, D. (2015) Advances in liquid chromatography–
9 high-resolution mass spectrometry for quantitative and qualitative environmental analysis.
10 *Analytical and bioanalytical chemistry* 407(21) 6289-6299.
- 11 [54] Hird, SJ.. Lau. BP-Y. Schuhmacher, R. Krska, R. (2014) Liquid chromatography-mass spectrometry
12 for the determination of chemical contaminants in food. *Trends in Analytical Chemistry*. 59; 59-72
- 13 [55] Zhang, Y. Li, Q. Lu, Y. Jones, K Sweetman, AJ. (2016). Hexabromocyclododecanes (HBCDDs) in
14 surface soils from coastal cities in North China: Correlation between diastereoisomer profiles and
15 industrial activities. *Chemosphere*. 148; 504-510.
- 16 [56] Managaki, S. Enomoto, I. Masunaga, S. (2012) Sources and distribution of hexabromocyclododecanes
17 (HBCDDs) in Japanese river sediment. *Journal of Environmental Monitoring* 14(3); 901-907
- 18 [57] Alder, L., Greulich, K., Kempe, G., Vieth, B., (2006). Residue analysis of 500 high priority pesticides:
19 better by GC–MS or LC–MS/MS? *Mass Spectrom. Rev*, 25, pp.838–865.
- 20 [58] Raina, R. (2011) Chemical analysis of pesticides using GC/MS, GC/MS/MS, and LC/MS/MS, in: M.
21 Stoycheva (Ed.), *Pesticides – Strategies for Pesticides Analysis*, InTech.
- 22 [59] Tienstra, M., Portolés, T., Hernández, F., Mol, JGJ., (2015). Fast gas chromatographic residue
23 analysis in animal feed using split injection and atmospheric pressure chemical ionisation tandem
24 mass spectrometry. *Journal of Chromatography A*. 1422; 289-298

25

- 1 [60] Covaci, A. Voorspoels, S. Ramos, L. Neels, H. Blust, R (2007) Recent developments in the analysis
2 of brominated flame retardants and brominated natural compounds. *J. Chromatogr. A.* 1153 (1-2);
3 145-171.
- 4 [61] Zacs, D., Rjabova, J., Pugajeva, I. Nakurte, I. Viksna, A., Bartkevics, V. (2014). Ultra high
5 performance liquid chromatography–time-of-flight high resolution mass spectrometry in the
6 analysis of hexabromocyclododecane diastereomers: Method development and comparative
7 evaluation versus ultra high performance liquid chromatography coupled to Orbitrap high
8 resolution mass spectrometry and triple quadrupole tandem mass spectrometry. *Journal of*
9 *Chromatography A.* 1366, 73-83
- 10 [62] Sales, C., Portoles, T., Sancho, JV., Abad, E., Abalos, M., Saulo, J., Fiedler, H., Gomara, B., Beltran,
11 J., 2016. Potential of gas chromatography-atmospheric pressure chemical ionization-tandem mass
12 spectrometry for screening and quantification of hexabromocyclododecane. *Analytical and*
13 *bioanalytical chemistry*, 408, 449-459.
- 14 [63] Kolic, TM., Shen, L., MacPherson, K., Fayez, L., Gobran, T., Helm, PA., Marvin, CH., Arsenault, G.,
15 Reiner, EJ., (2009). The Analysis of Halogenated Flame Retardants by GC–HRMS in
16 Environmental Samples. *Journal of Chromatographic Science*, 47, 83-91
- 17 [64] Jahnke, A. Berger, U. (2009) Trace analysis of per- and polyfluorinated alkyl substances in various
18 matrices-How do current methods perform? *J. Chromatogr. A.* 1216 (3); 410-421
- 19 [65] Chu, S., Letcher, RJ., McGoldrick, DJ., Backus, SM., (2016). A New Fluorinated Surfactant
20 Contaminant in Biota: Perfluorobutane Sulfonamide in Several Fish Species. *Environmental*
21 *Science and Technology*, 50, 669–675
- 22 [66] D'Agostino, LA., Mabury, SA., 2014. Identification of Novel Fluorinated Surfactants in Aqueous Film
23 Forming Foams and Commercial Surfactant Concentrates. *Environmental Science and*
24 *Technology*. 48(1), 121-129
- 25 [67] Guo, R., Reiner, EJ., Bhavsar, SP., Helm, PA., Mabury, SA., Braekevelt, E., Tittlemier, SA. (2012).
26 Determination of polyfluoroalkyl phosphoric acid diesters, perfluoroalkyl phosphonic acids,

- 1 perfluoroalkyl phosphinic acids, perfluoroalkyl carboxylic acids, and perfluoroalkane sulfonic
2 acids in lake trout from the Great Lakes region, *Analytical and Bioanalytical Chemistry*, 404(9),
3 2699-2709
- 4 [68] Lui, Y., Pereira, ADS., Martin, JW., (2015). Discovery of C5–C17 Poly- and Perfluoroalkyl
5 Substances in Water by In-Line SPE-HPLC-Orbitrap with In-Source Fragmentation Flagging.
6 *Analytical Chemistry*, 87(8) 4260-4268
- 7 [69] Phillips, JB. Lui, Z. (1991) *Comprehensive Two-Dimensional Gas Chromatography using an On-*
8 *Column Thermal Modulator Interface. J. Chromatogr. Sci.* 29 (6): 227-231
- 9 [70] Adahchour, M. Beens, J. Brinkman, UAT. (2008) Recent developments in the application of
10 comprehensive two-dimensional gas chromatography. *J. Chromatogr. A.* 1186: 67-108
- 11 [71] Bordajandi LR, Ramos JJ, Sanz J, Gonzalez MJ, Ramos L, (2008) Comprehensive two-dimensional
12 gas chromatography in the screening of persistent organohalogenated pollutants in environmental
13 samples. *J. Chromatogr. A.* 1186: 312-324.
- 14 [72] Mondello, L. Tranchida, PQ. Dugo, P. Dugo, G. (2008) Comprehensive two-dimensional gas
15 chromatography-mass spectrometry: a review, *Mass Spectrom. Rev.* 27: 101-124.
- 16 [73] Murray, JA. (2012) Qualitative and quantitative approaches in comprehensive two-dimensional gas
17 chromatography. *J. Chromatogr. A.* 1261, pp. 58-68.
- 18 [74] Meinert, C. Meierhenrich, UJ. (2012) A new dimension in separation science: Comprehensive two-
19 dimensional gas chromatography, *Angewandte Chemie - International Edition*, 51(42): 10460-
20 10470.
- 21 [75] Mostafa, A. Edwards, M. Górecki, T. (2012) Optimization aspects of comprehensive two-dimensional
22 gas chromatography. *J. Chromatogr. A.* 1255: 38-55.
- 23 [76] Marriott, PJ. Chin, S-T. Maikhunthod, B. Schmarr, H-G. Bieri, S. (2012) Multidimensional gas
24 chromatography. *Trends Anal. Chem.* 34: 1-20
- 25 [77] Ramos L (2009) *Comprehensive Two Dimensional Gas Chromatography.* Elsevier, Oxford, UK

- 1 [78] Seeley, JV. Seeley, SK. (2013) Multidimensional gas chromatography: Fundamental advances and
2 new applications, *Anal. Chem.* 85: 557-578.
- 3 [79] Danielsson, C. Wiberg, K. Korytar, P. Bergek, S. Brinkman, UATh. Haglund, P. (2005) Trace analysis
4 of polychlorinated dibenzo-p-dioxins, dibenzofurans and WHO polychlorinated biphenyls in food
5 using comprehensive two-dimensional gas chromatography with electron-capture detection. *J.*
6 *Chromatogr. A.* 1086: 61-70
- 7 [80] Focant, JF. Reiner, EJ. MacPherson, K. Kolic, T. Sjodin, A. Patterson Jr, DG. Reese, SL. Dorman, FL.
8 Cochran, J. (2004a), Measurement of PCDDs, PCDFs, and non-ortho-PCBs by comprehensive
9 two-dimensional gas chromatography-isotope dilution time-of-flight mass spectrometry (GC x GC-
10 IDTOFMS). *Talanta.* 63: 1231-1240
- 11 [81] Focant, JF. Sjodin, A. Patterson Jr, DG. (2004b), Improved separation of the 209 polychlorinated
12 biphenyl congeners using comprehensive two-dimensional gas chromatography-time-of-flight
13 mass spectrometry. *J. Chromatogr. A.* 1040: 227-23
- 14 [82] Focant, JF. Sjodin, A. Turner, WE. Patterson Jr, DG. (2004c), Measurement of selected
15 polybrominated diphenyl ethers, polybrominated and polychlorinated biphenyls, and
16 organochlorine pesticides in human serum and milk using comprehensive two-dimensional gas
17 chromatography isotope dilution time-of-flight mass spectrometry. *Anal. Chem.* 76: 6313-6320
- 18 [83] Haglund, P. Korytar, P. Danielsson, C. Diaz, J. Wiberg, K. Leonards, P. Brinkman, IAT. de Boer, J.
19 (2008) GCxGC-ECD: a promising method for the determination of dioxins and dioxin-like PCBS
20 in food and feed. *Anal. Bioanal. Chem.* 390: 1815-1827.
- 21 [84] Korytar, P. Danielsson, C. Leonards, PEG. Haglund, P. de Boer, J. Brinkman, UATh. (2004)
22 Separation of seventeen 2,3,7,8-substituted polychlorinated dibenzo-p-dioxins and dibenzofurans
23 and 12 dioxin-like polychlorinated biphenyls by comprehensive two-dimensional gas
24 chromatography with electron-capture detection. *J.Chromatogr. A* 1038: 189-199.
- 25 [85] Korytar, P. Parera, J. Leonards, PEG. Santos, FJ. de Boer, J. Brinkman, UAT. (2005) Characterization
26 of polychlorinated n-alkanes using comprehensive two-dimensional gas chromatography-electron-

- 1 capture negative ionisation time-of-flight mass spectrometry. *Journal of Chromatography A* 1086
2 (1-2):71-82
- 3 [86] Muscalu, AM. Reiner, EJ. Liss, SN. Chen, T. Ladwig, G. Morse, D. (2011) A routine accredited
4 method for the analysis of polychlorinated biphenyls, organochlorine pesticides, chlorobenzenes
5 and screening of other halogenated organics in soil, sediment and sludge by GCxGC- μ ECD.
6 *Anal. Bioanal. Chem.* 401 (8):2403-2413.
- 7 [87] Megson, D. Brown, TA. Johnson, GW. O'Sullivan. G. Bicknell, AWJ. Votier, SC. Lohan, MC.
8 Comber, S. Kalin, R. Worsfold, PJ. (2014) Identifying the provenance of Leach's storm petrels in
9 the North Atlantic using polychlorinated biphenyl signatures derived from comprehensive two-
10 dimensional gas chromatography with time-of-flight mass spectrometry. *Chemosphere.* 114: 195-
11 202
- 12 [88] Sampat, A., Lopatka, M., Sjerps, M., Vivo-Truyols, G., Schoenmakers, P., van Asten, A., (2016).
13 Forensic potential of comprehensive two-dimensional gas chromatography. *Trends in Analytical*
14 *Chemistry.* 80, 345–363
- 15 [89] Rimayi, C., Chimuka, L., Odusanya, D., de Boer, J., Weiss, J., 2016. Distribution of 2,3,7,8-substituted
16 polychlorinated dibenzo-p-dioxin and polychlorinated dibenzofurans in the Jukskei and Klip/Vaal
17 catchment areas in South Africa. *Chemosphere.* 145, 314-321.
- 18 [90] Guo Y, Kannan K, (2015) Analytical Methods for the Measurement of Legacy and Emerging Persistent
19 Organic Pollutants in Complex Sample Matrices, In Zeng EY, *Comprehensive Analytical*
20 *Chemistry, Persistent Organic Pollutants (POPs): Analytical Techniques, Environmental Fate and*
21 *Biological Effects, Elsevier Vol 67, 1-56*
- 22 [91] M. Petrovic, M. Farre, M.L. de Alda, S. Perez, C. Postigo, M. Kock, J. Radjenovic, M. Gros, D.
23 Barcelo, (2010) Recent trends in the liquid chromatography-mass spectrometry analysis of organic
24 contaminants in environmental samples, *J. Chromatogr. A.* 1217 4004–4017.
- 25 [92] F. Hernández, J.V. Sancho, M. Ibáñez, E. Abad, T. Portolés, L. Mattioli, Current use of high-resolution
26 mass spectrometry in the environmental sciences, (2012) *Anal. Bioanal. Chem.* 403. 1251–1264.

- 1 [93] M. Farré, L. Kantiani, M. Petrovic, S. Pérez, D. Barceló, (2012) Achievements and future trends in the
2 analysis of emerging organic contaminants in environmental samples by mass spectrometry and
3 bioanalytical techniques, *J. Chromatogr. A* 1259 86–99.
- 4 [94] Harrison, AG. (1992). *Chemical Ionization Mass Spectrometry*, (Second Edition). CRC Press
- 5 [95] Budde, WL. (in press). Ron Hites: Gas Chromatography/Mass Spectrometry Pioneer and the Great
6 Detector Debate. *Environ. Sci. Technol.* (DOI: 10.1021/es505993z)
- 7 [96] Zimmermann, R. Welthagen, W. Groger, T. (2008) Photo-ionisation mass spectrometry as detection
8 method for gas chromatography. Optical selectivity and multidimensional comprehensive
9 separations. *Journal of Chromatography A*. 1184; 296–308
- 10 [97] Li, DX. Gan, L. Bronja, A. Schmitz, OJ. (2015) Gas chromatography coupled to atmospheric pressure
11 ionization mass spectrometry (GC-API-MS): Review. *Anal. Chim. Acta*. 891, 43-61
- 12 [98] van Bavel, B. Geng, D. Cherta, L. Náchter-Mestre, J. Portolés, T. Ábalos, M. Sauló, J. Abad, E.
13 Dunstan, J. Jones, R. Kotz, A. Winterhalter, H. Malisch, R. Traag, W. Hagberg, J. Ericson Jogsten,
14 I. Beltran, J. and Hernández, F. (2015) Atmospheric-Pressure Chemical Ionization Tandem Mass
15 Spectrometry (APGC/MS/MS) an Alternative to High-Resolution Mass Spectrometry
16 (HRGC/HRMS) for the Determination of Dioxins. *Anal. Chem.* 87: 9047–9053
- 17 [99] Organtini, KL. Haimovici, L. Jobst, K. Reiner, EJ. Ladak, A. Stevens, D. Cochran, JW. Dorman, FL.
18 (2015) Comparison of Atmospheric Pressure Ionization Gas Chromatography-Triple Quadrupole
19 Mass Spectrometry to Traditional High-Resolution Mass Spectrometry for the Identification and
20 Quantification of Halogenated Dioxins and Furans. *Anal. Chem.* 87: 7902-7908.
- 21 [100] Ballesteros-Gómez, A. de Boer, J. Leonards, PEG. (2013) Novel Analytical Methods for Flame
22 Retardants and Plasticizers Based on Gas Chromatography, Comprehensive Two-Dimensional Gas
23 Chromatography, and Direct Probe Coupled to Atmospheric Pressure Chemical Ionization-High
24 Resolution Time-of-Flight-Mass Spectrometry. *Anal. Chem.* 85: 9572-9580
- 25 [101] Portoles, T. Sales, C. Gómara, B. Sancho, JV. Beltrán, J. Herrero, L. González, MJ. Hernández, F.
26 (2015) Novel Analytical Approach for Brominated Flame Retardants Based on the Use of Gas

- 1 Chromatography-Atmospheric Pressure Chemical Ionization-Tandem Mass Spectrometry with
2 Emphasis in Highly Brominated Congeners. *Anal. Chem.* 87: 9892-9899.
- 3 [102] Lobodin, VV. Rodgers, RP. Marshal, AG. (2012) Chapter 18 Petroleomics and the Analysis of
4 Complex Organic Mixtures with Fourier Transform Ion Cyclotron Resonance. in *Comprehensive*
5 *Environmental Mass Spectrometry (Advanced Topics in Environmental Science): 4* (Lebedev, A.T.
6 ed). ILM Publications.
- 7 [103] Korcz, W. Strucinski, P. Goralczyk, K. Hernik, A. Lyczewska, M. Czaja, K. Matuszak, M.
8 Minorczyk, M. Ludwicki, J.K. 2014. Development and validation of a method for determination of
9 selected polybrominated diphenyl ether congeners in household dust. *Roczniki Panstwowego*
10 *Zakladu Higieny.* 65,2. 93-100
- 11 [104] Qi, H. Li, W.L. Zhu, N.Z. Ma, W.L. Liu, L.Y. Zhang, F. Li, Y.F. (2014) Concentrations and sources
12 of polycyclic aromatic hydrocarbons in indoor dust in Chinal. *Science of the Total Environment.*
13 491. 100-107
- 14 [105] Kalachova, K. Cajka, T. Sandy, C. Hajslova, J. Pulkrabova, J. 2013. High throughput sample
15 preparation in combination with gas chromatography coupled to triple quadrupole tandem mass
16 spectrometry (GC-MS/MS): A smart procedure for (ultra)trace analysis of brominated flame
17 retardants in fish. *Talanta.* 105. 109-116
- 18 [106] Brandsma, SH., Leonards, PEG., Leslie, HA., de Boer, J., (2015). Tracing organophosphorus and
19 brominated flame retardants and plasticizers in an estuarine food web. *Science of the Total*
20 *Environment,* 505, 22–31
- 21 [107] Myers, A.L., Watson-Leung, T., Jobst, K.J., Shen, L., Besevic, S., Organtini, K., Dorman, F.L.,
22 Mabury, S.A. and Reiner, E.J., 2014. Complementary nontargeted and targeted mass spectrometry
23 techniques to determine bioaccumulation of halogenated contaminants in freshwater species.
24 *Environmental science & technology,* 48(23), pp.13844-13854.
- 25 [108] Patterson Jr, DG. Welch, SM. Turner, WE. Sjodin, A. Focant, JF. (2011) Cryogenic zone compression
26 for the measurement of dioxins in human serum by isotope dilution at the attogram level using

- 1 modulated gas chromatography coupled to high resolution magnetic sector mass spectrometry. J.
2 Chromatogr. A. 1218, 3274-3281.
- 3 [109] Ortiz- Almirall X, Pena-Abaurrea M, Jobst KJ, Reiner EJ, (2016) Non-Targeted Analysis of Persistent
4 Organic Pollutants by Mass Spectrometry and GC×GC, In Perez S, Eichhorn P, Barcelo D, (eds),
5 Comprehensive Analytical Chemistry Applications of Time-of-Flight and Orbitrap Mass
6 Spectrometry in Environmental, Food, Doping, and Forensic, Elsevier, Vol 71 (in Press)
- 7 [110] Fernando, S.; Jobst, K. J.; Taguchi, V. Y.; Helm, P. A.; Reiner, E. J.; McCarry, B. E. (2014)
8 Identification of the halogenated compounds resulting from the 1997 Plastimet Inc. fire in
9 Hamilton, Ontario, using comprehensive two-dimensional gas chromatography and (ultra)high
10 resolution mass spectrometry. Environmental Science and Technology. 48, 10656– 10663.
- 11 [111] Peterson, AC. McAlister, GC. Quarmby, GC. Griep-Raming, J. Coon, JJ. (2010). Development and
12 Characterization of a GC-Enabled QLT-Orbitrap for High-Resolution and High-Mass Accuracy
13 GC/MS. Analytical Chemistry. 82, 8618-8628.
- 14 [112] Zhao, XS., Zhou, YK., Kong, WJ., Gong, B., Chen, DL., Wei, JH., Yang, MH. (2016). Multi-residue
15 analysis of 26 organochlorine pesticides in *Alpinia oxyphylla* by GC-ECD after solid phase
16 extraction and acid cleanup. Journal of Chromatography B, 1017, 211-220
- 17 [113] Chawla, S., Patel, HK., Vaghela, KM., Pathan, FK., Gor, HN., Patel, AR., Shah, PG., (2016).
18 Development and validation of multiresidue analytical method in cotton and groundnut oil for 87
19 pesticides using low temperature and dispersive cleanup on gas chromatography and liquid
20 chromatography-tandem mass spectrometry. Analytical and Bioanalytical Chemistry, 408(3), 983-
21 997
- 22 [114] Tette, PAS., Guidi, LR., Gloria, MBdA., Fernandes, C. (2016). Pesticides in honey: A review on
23 chromatographic analytical methods. Talanta. 149, 124-141.
- 24 [115] Shendy, AH., Al-Ghobashy, MA., Mohammed, MN., Alla, SAG., Lofty, HM. 2016. Simultaneous
25 determination of 200 pesticide residues in honey using gas chromatography–tandem mass

- 1 spectrometry in conjunction with streamlined quantification approach. *Journal of Chromatography*
2 A, 1427, 142-160
- 3 [116] Wong, JW. Zhang, K., Tech, K., Hayward, DG., Makovi, CM., Krynitsky, AJ., Schenck, FJ.,
4 Banerjee, K., Dasgupta, S., Brown, D. 2010. Multiresidue Pesticide Analysis in Fresh Produce by
5 Capillary Gas Chromatography-Mass Spectrometry/Selective Ion Monitoring (GC-MS/SIM) and-
6 Tandem Mass Spectrometry (GC-MS/MS). *Journal of Agricultural and Food Chemistry*, 58(10),
7 5868-5883
- 8 [117] Gomez-Ramos, MD., Rajski, L., Heinzen, H., Fernandez-Alba, AR. 2015. Liquid chromatography
9 Orbitrap mass spectrometry with simultaneous full scan and tandem MS/MS for highly selective
10 pesticide residue analysis. *Analytical and Bioanalytical Chemistry*, 407(21), 6317-6326
- 11 [118] Kim, MS. Kang, TW. Pyo, H. Yoon, J. Choi, K., Hong, J. (2008). Determination of organochlorine
12 pesticides in sediment using graphitized carbon black solid-phase extraction and gas
13 chromatography/mass spectrometry. *Journal of Chromatography A*, 1208(1-2), 25-33
- 14 [119] Megson, D., Ortiz, X., Jobst, KJ., Reiner, EJ., Mulder, MFA., Balouet, JC. (2016). A comparison of
15 fresh and used aircraft oil for the identification of toxic substances linked to aerotoxic syndrome.
16 *Chemosphere*. 158, 116-123
- 17 [120] Silva Barni, MF. Ondarza, PM. Gonzalez, M. Da Cuna, R. Meijide, F. Grosman, F. Sanzano, P. Lo
18 Nostro, FL. Miglioranza, KSB. (2016). Persistent organic pollutants (POPs) in fish with different
19 feeding habits inhabiting a shallow lake ecosystem. *Science of the Total Environment*, 550, 900–
20 909
- 21 [121] Ben Salem, F., Ben Said, O., Duran, R., Monperrus, M. (2016) Validation of an Adapted QuEChERS
22 Method for the Simultaneous Analysis of Polycyclic Aromatic Hydrocarbons, Polychlorinated
23 Biphenyls and Organochlorine Pesticides in Sediment by Gas Chromatography–Mass
24 Spectrometry. *Bull Environ Contam Toxicol*. 96(5); 678-84.
- 25 [122] Ochiai, N. Ieda, T. Sasamoto, K., Takazawa, Y. Hashimoto, S. Fushimi, A. Tanabe, K. (2011). Stir
26 bar sorptive extraction and comprehensive two-dimensional gas chromatography coupled to high-

- 1 resolution time-of-flight mass spectrometry for ultra-trace analysis of organochlorine pesticides in
2 river water. *Journal of Chromatography A*, 1218(39), 6851-6860
- 3 [123] Salihovic, S. Mattioli, L. Lindstrom, G. Lind, L. Lind, PM. van Bavel, B. (2012). A rapid method for
4 screening of the Stockholm Convention POPs in small amounts of human plasma using SPE and
5 HRGC/HRMS. *Chemosphere*, 86(7); 747-753
- 6 [124] Farre, M. Pico, Y. Barcelo, D. (2014). Application of ultra-high pressure liquid chromatography linear
7 ion-trap orbitrap to qualitative and quantitative assessment of pesticide residues. *Journal of*
8 *Chromatography A*, 1328, 66-79
- 9 [125] Castells, P., Parera, J., Santos, FJ., Galceran, MT. (2008). Occurrence of polychlorinated
10 naphthalenes, polychlorinated biphenyls and short-chain chlorinated paraffins in marine sediments
11 from Barcelona (Spain). *Chemosphere*, 70(9); 1552-1562
- 12 [126] Zaater, M. Tahboub, Y. Qasrawy, S. (2005). Monitoring of polychlorinated biphenyls in surface water
13 using liquid extraction, GC/MS, and GC/ECD. *Analytical Letters*. 38(13), 2231-2245.
- 14 [127] Montory, M. Habit, E. Fernandez, P. Gribmalt, JO. Barra, R. (2010). PCBs and PBDEs in wild
15 Chinook salmon (*Oncorhynchus tshawytscha*) in the Northern Patagonia, Chile. *Chemosphere*, 78,
16 1193–1199
- 17 [128] Barco-Bonilla, N. Nieto-Garcia, AJ. Romero-Gonzalez, R. Vidal, JLM. Frenich, AG. (2015).
18 Simultaneous and highly sensitive determination of PCBs and PBDEs in environmental water and
19 sediments by gas chromatography coupled to high resolution magnetic sector mass spectrometry.
20 *Analytical Methods*, 7(7), pp.3036-3047
- 21 [129] Naccha, L. Alanis, G. Torres, A. Abad, E. Abalos, M. Rivera, J. Heyer, L. Morales, A. Waksman, N.
22 (2010). Dioxins in beef samples from Mexico using a low resolution GC/MS screening method,
23 *Food Additives & Contaminants: Part B*, 3(1), 64-72
- 24 [130] de Vos, J. Gorst-Allman, P. Rohwer, E. (2011). Establishing an alternative method for the quantitative
25 analysis of polychlorinated dibenzo-p-dioxins and polychlorinated dibenzofurans by

- 1 comprehensive two dimensional gas chromatography–time-of-flight mass spectrometry for
2 developing countries. *Journal of Chromatography A*, 1218, 3282–3290
- 3 [131] Haimovici L, Reiner EJ, Besevic S, Jobst KJ, Robson M, Kolic T, MacPherson K, (2016) A Modified
4 QuEChERS Approach for the Screening of Dioxins and Furans in Sediments, *Anal Bioanal Chem*,
5 408(15), pp.4043-4054
- 6 [132] Hu, XZ., Hu, DC., Chen, W., Wu, B., Lin, CJ., 2015. Simultaneous Determination of Methoxylated
7 Polybrominated Diphenyl Ethers and Polybrominated Diphenyl Ethers in Water, Soil and Sediment
8 from China by GC-MS. *Journal of Chromatographic Science*, 53(8), 1239-1249
- 9 [133] Kalachova, K., Pulkrabova, J, Cajka, T, Drabova, L, Hajslova, J, (2012). Implementation of
10 comprehensive two-dimensional gas chromatography-time-of-flight mass spectrometry for the
11 simultaneous determination of halogenated contaminants and polycyclic aromatic hydrocarbons in
12 fish. *Analytical and Bioanalytical Chemistry*, 403(10), 2813-2824
- 13 [134] Song, SJ. Dai, XH. Wang, WH. He, YJ. Lui, Z. Shao, MW. (2015). Optimization of selective
14 pressurized liquid extraction and ultrasonication-assisted QuEChERS methods for the
15 determination of polybrominated diphenyl ethers in sediments. *Analytical Methods*. 7(22), 9542-
16 9548
- 17 [135] Fromme, H. Hilher, B. Kopp, E. Miserok, M. Volkel, W. (2014). Polybrominated diphenyl ethers
18 (PBDEs), hexabromocyclododecane (HBCD) and “novel” brominated flame retardants in house
19 dust in Germany. *Environment International*, 64, 61–68
- 20 [136] Lankova, D. Svarcova, A. Kalachova, K. Lacina, O. Pulkrabova, J. (2015) Multi-analyte method for
21 the analysis of various organohalogen compounds in house dust. *Analytica Chimica Acta*, 854, 61–
22 69
- 23 [137] Zacs, D. Bartkevics, V. (2015). Analytical capabilities of high performance liquid chromatography e
24 Atmospheric pressure photoionization e Orbitrap mass spectrometry (HPLC-APPI-Orbitrap-MS)
25 for the trace determination of novel and emerging flame retardants in fish. *Analytica Chimica Acta*.
26 898, 60-72.

- 1 [138] Darnerud, PO. Aune, M. Larsson, L. Lignell, S. Mutshatshi, T. Okonkwo, J. Botha, B. Agyei, N.
2 (2011). Levels of brominated flame retardants and other persistent organic pollutants in breast milk
3 samples from Limpopo province, South Africa. *Science of the Total Environment*, 409, 4048–4053
- 4 [139] Chokwe, TB. Okonkwo, JO. Sibali, L. Ncube, EJ. (2012). Optimization and Simultaneous
5 Determination of Alkyl Phenol Ethoxylates and Brominated Flame Retardants in Water after SPE
6 and Heptafluorobutyric Anhydride Derivatization followed by GC/MS. *Chromatographia*, 75,
7 1165–1176
- 8 [140] Abdallah, MAE. Ibarra, C., Neels, H., Harrad, S., Covaci, A. (2008). Comparative evaluation of liquid
9 chromatography–mass spectrometry versus gas chromatography–mass spectrometry for the
10 determination of hexabromocyclododecanes and their degradation products in indoor dust. *Journal*
11 *of Chromatography A*, 1190, 333–341
- 12 [141] Sun, C. Sun, H. Lai, Y. Zhang, J. Cai, Z. (2011), Liquid Chromatography/Mass Spectrometry Method
13 for Determination of Perfluorooctane Sulfonyl Fluoride upon Derivatization with Benzylamine.
14 *Analytical Chemistry*, 83, 5822–5826
- 15 [142] Padilla-Sanchez, JA. Haug, JS. (2016). A fast and sensitive method for the simultaneous analysis of
16 a wide range of per- and polyfluoroalkyl substances in indoor dust using on-line solid phase
17 extraction-ultrahigh performance liquid chromatography-time-of-flight-mass spectrometry.
18 *Journal of Chromatography A*, 1445, 36-45.
- 19 [143] Campo, J. Lorenzo, M. Perez, F. Pico, Y. Farre, M. Barcelo, D. (2016). Analysis of the presence of
20 perfluoroalkyl substances in water, sediment and biota of the Jucar River (E Spain). Sources,
21 partitioning and relationships with water physical characteristics. *Environmental Research*, 147,
22 503-512.
- 23 [144] Munoz, G. Vo Duy, S. Budzinski, H. Labadie, P. Lui, J. Sauve, S. (2015). Quantitative analysis of
24 poly- and perfluoroalkyl compounds in water matrices using high resolution mass spectrometry:
25 Optimization for a laser diode thermal desorption method. *Analytica Chimica Acta* 881, 98–106.

- 1 [145] Johnson, GW. Quensen III, JF., Chiarenzelli, JR. Coreen Hamilton, M. (2006). Polychlorinated
2 Biphenyls. In: Morrison, R.D. & Murphy, B.L. (eds.) Environmental Forensics Contaminant
3 Specific Guide. Academic Press.
- 4 [146] Balouet, JC. Smith, KT. Vroblesky, D. Oudijk, G. (2009) Use of Dendrochronology and
5 Dendrochemistry in Environmental Forensics: Does It Meet the Daubert Criteria? Environmental
6 Forensics. 10: 268–276
- 7 [147] Brillis, GM. Worthington, JC. and Wait, AD. (2000) Quality science in the courtroom: U.S. EPA data
8 quality and peer review policies and procedures compared to the Daubert factors. Environmental
9 Forensics. 1: 197–203.
- 10 [148] Kanner, A. (2007). Meeting the burden for admissibility of environmental forensic testimony.
11 Environmental Forensics. 8: 19–23.
- 12 [149] Ries, DG. Burns, RL. (2005) Expert opinions in environmental litigation: Gatekeeping twelve years
13 after Daubert. Environmental Forensics. 5(3): 219–229.
- 14 [150] HimaBindu LMR, Parameswari AS, Gopinath C, (2013), A review on GC-MS method development
15 and validation, International Journal of Pharmaceutical Quality Assurance 4;3. 42-51
- 16 [151] Frame, GM. (2001) The Current State-of-the-Art of Comprehensive, Quantitative, Congener-Specific
17 PCB Analysis, and What We Now Know about the Distributions of Individual Congeners in
18 Commercial Aroclor Mixtures. In: Robertson LW, Hansen LG, editors. PCBs Recent Advances in
19 Environmental Toxicology and Health Effects. The University Press of Kentucky, Lexington, KY
- 20 [152] Creaser, CS. Wood, MD. Alcock, R. Copplestone, D. Crook, PJ. Barraclough, D. (2007)
21 Environmental concentrations of polychlorinated dibenzo-p-dioxins and polychlorinated
22 dibenzofurans in UK soil and herbage, UKSHS Report No. 10, Bristol, Environment Agency.
- 23 [153] Balouet, J-C. Gallion, F. Martelain, J, Megson, D. O’Sullivan, G. 2015. Forensic investigations in
24 complex pollution cases involving pcbs, dioxins and furans: potential pitfalls and tips. In
25 O’Sullivan, G and Megson, D. (eds) Environmental Forensics, Proceedings of the 2014 conference.
26 RSC publishing. Cambridge, UK.

- 1 [154] Richman L, Haimovici L, Kolic T, Besevic S, Reiner E, (2016) Monitoring Re-Suspension and
2 Transport of Dioxin Contaminated Sediment to Evaluate the Recovery of a Shallow Urban Creek
3 Post Sediment Remediation, *Journal of Environmental Protection*, 7, 453-466
- 4 [155] Bertazzi, P.A. Bernucci, I; Brambilla, G. Consonni, D. Pesatori, AC. 1996. The Seveso studies on
5 early and long-term effects of dioxin exposure. *Environmental Health Perspectives*. 106;2. 625-633
- 6 [156] Brown, P 1992. Popular epidemiology and toxic-waste contamination - lay and professional ways of
7 knowing. *Journal of health and social behavior*. 33;3. 267-281
- 8 [157] Gago-Ferrero, P., Schymanski, E.L., Bletsou, A.A., Aalizadeh, R., Hollender, J. and Thomaidis, N.S.,
9 (2015). Extended suspect and non-target strategies to characterize emerging polar organic
10 contaminants in raw wastewater with LC-HRMS/MS. *Environmental science & technology*,
11 49(20).12333-12341.
- 12 [158] Schymanski, E.L., Singer, H.P., Slobodnik, J., Ipolyi, I.M., Oswald, P., Krauss, M., Schulze, T.,
13 Haglund, P., Letzel, T., Grosse, S. and Thomaidis, N.S., (2015) Non-target screening with high-
14 resolution mass spectrometry: critical review using a collaborative trial on water analysis.
15 *Analytical and bioanalytical chemistry*, 407(21), pp.6237-6255.
- 16 [159] Strynar, M., Dagnino, S., McMahan, R., Liang, S., Lindstrom, A., Andersen, E., McMillan, L.,
17 Thurman, M., Ferrer, I. and Ball, C., (2015) Identification of Novel Perfluoroalkyl Ether Carboxylic
18 Acids (PFECAs) and Sulfonic Acids (PFESAs) in Natural Waters Using Accurate Mass Time-of-
19 Flight Mass Spectrometry (TOFMS). *Environmental science & technology*, 49(19), pp.11622-
20 11630.
- 21 [160] Bletsou, A.A., Jeon, J., Hollender, J., Archontaki, E. and Thomaidis, N.S., (2015) Targeted and non-
22 targeted liquid chromatography-mass spectrometric workflows for identification of transformation
23 products of emerging pollutants in the aquatic environment. *TrAC Trends in Analytical Chemistry*,
24 66, pp.32-44.
- 25 [161] Parette, R. and Pearson, W.N. 2014. 2,4,6,8-Tetrachlorodibenzothiophene in the Newark Bay
26 Estuary: The likely source and reaction pathways. *Chemosphere*. 111. 157–163

- 1 [162] Sverko, E. Reiner, E.J. Tomy, G.T. McCrindle, R. Shen, L. Arsenault, G. Zaruk, D. MacPherson,
2 K.A. Marvin, C.H. Helm, P.A. McCarry, B.E. (2010) Compounds Structurally Related to
3 Dechlorane Plus in Sediment and Biota from Lake Ontario (Canada)
- 4 [163] Kendrick, E. (1963) A Mass Scale Based on $CH_2 = 14.0000$ for High Resolution Mass Spectrometry
5 of Organic Compounds, *Anal. Chem.* 35 (13) 2146-2154.
- 6 [164] Jobst, K.J. Shen, L. Reiner, E.J. Taguchi, V.Y. Helm, P.A. McCrindle, R. and Backus, S. (2013). The
7 use of mass defect plots for the identification of (novel) halogenated contaminants in the
8 environment. *Analytical and bioanalytical chemistry*, 405(10), 3289-3297.
- 9 [165] Pleil JD, Issacs KK, (2016) High-resolution mass spectrometry: basic principles for using exact mass
10 and mass defect for discovery analysis of organic molecules in blood, breath, urine and
11 environmental media, *J. Breath Res.* 10 012001
- 12 [166] Ubukata, M. Jobst, K.J. Reiner, E.J. Reichenbach, S.E. Tao, Q. Hang, J. Wu, Z. Dane, A.J. Cody,
13 R.B. (2015). Non-targeted analysis of electronics waste by comprehensive two-dimensional gas
14 chromatography combined with high-resolution mass spectrometry: Using accurate mass
15 information and mass defect analysis to explore the data. *Journal of Chromatography A*, 1395,152-
16 159.
- 17 [167] Duncan, K.D. Letourneau, D.R. Vandergrift, G.W. Jobst, K. Reiner, E. Gill, C.G. and Krogh, E.T.
18 (2016). A semi-quantitative approach for the rapid screening and mass profiling of naphthenic acids
19 directly in contaminated aqueous samples. *Journal of Mass Spectrometry*, 51(1), pp.44-52.
- 20 [168] Jobst, K.J. Pena-Abaurrea, M. Reiner, E.J. (2014), Strategies and Techniques for Identifying
21 Unknown Compounds in Environmental Samples. In *Environmental Forensics: Proceedings of the*
22 *2013 INEF Conference*. Royal Society of Chemistry.
- 23 [169] Myers, A.L. Jobst, K.J. Mabury, S.A. and Reiner, E.J. (2014). Using mass defect plots as a discovery
24 tool to identify novel fluoropolymer thermal decomposition products. *Journal of Mass*
25 *Spectrometry*, 49(4), pp.291-296.

26

1

2